



# MX from an Industrial Perspective


Dr. Ailsa Powell

Senior Industrial Liaison Scientist

[ailsa.powell@diamond.ac.uk](mailto:ailsa.powell@diamond.ac.uk)



# Types of industry and advantages of using MX

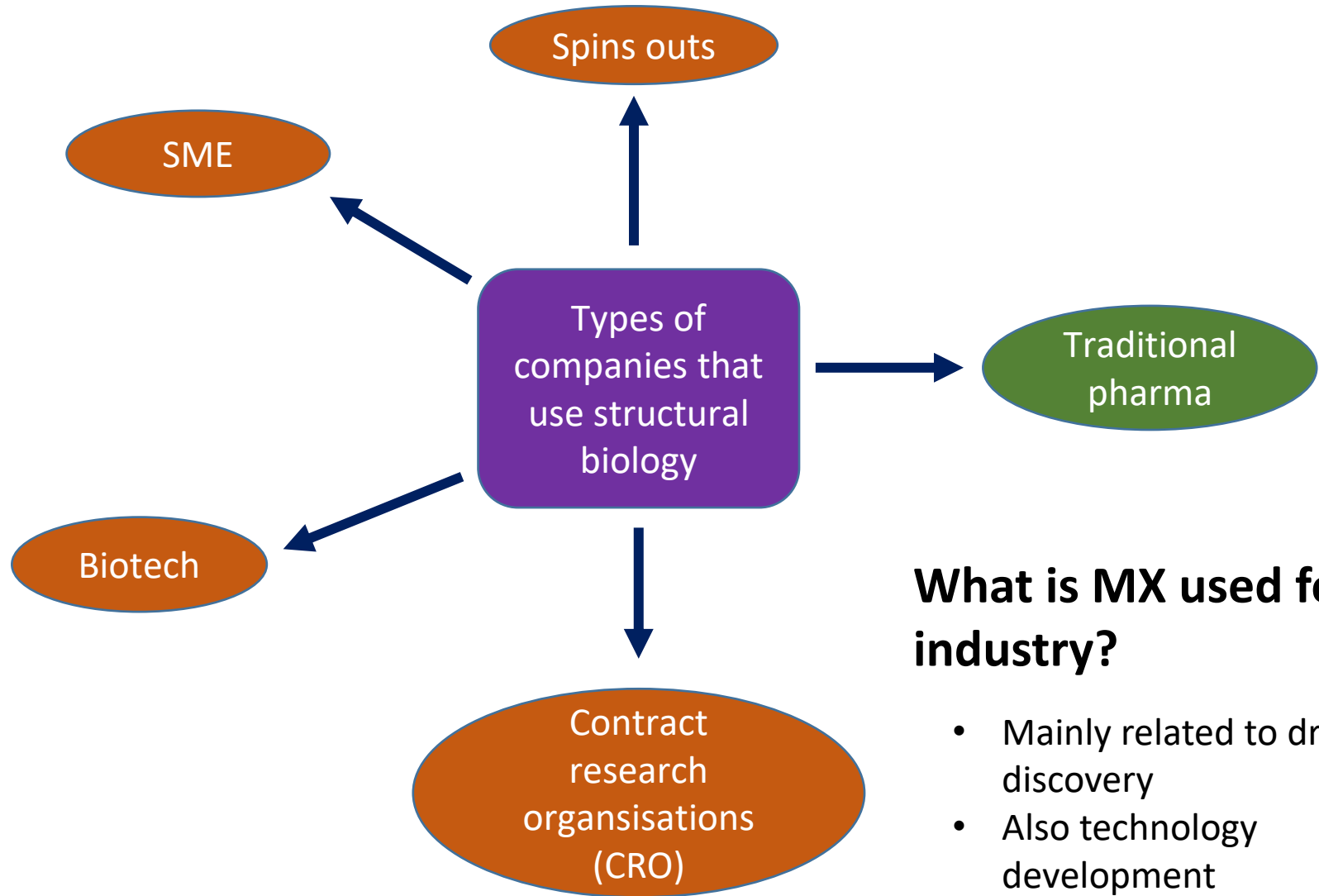


## What is MX used for in industry?


- Mainly related to drug discovery
- Also technology development

# Why use structural biology

Invaluable insights at the drug discovery phase

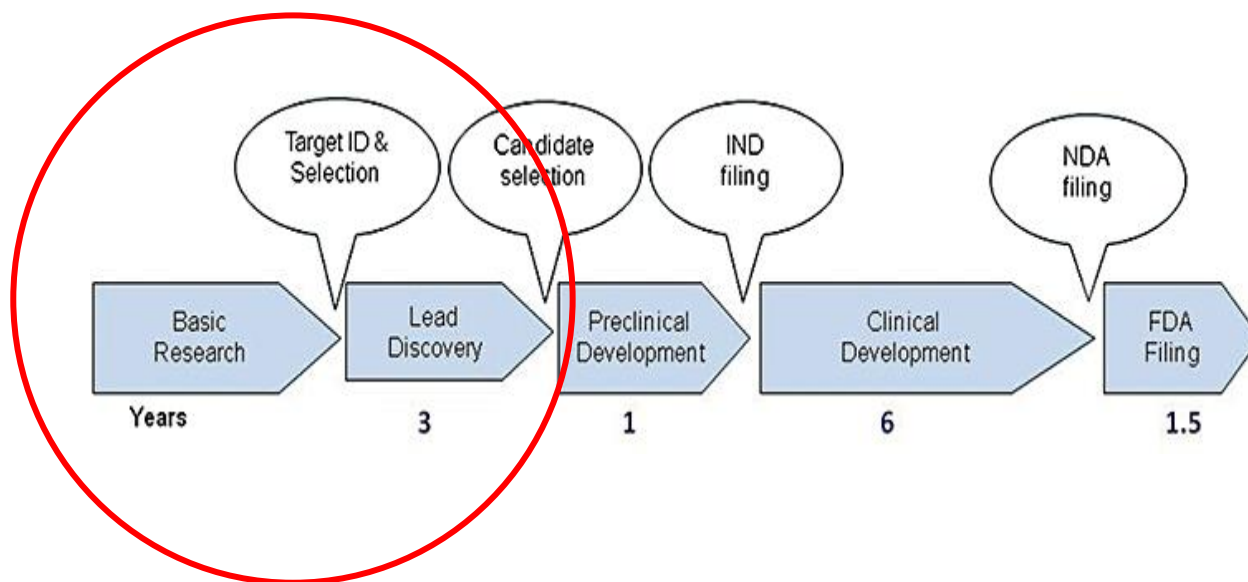

- Protein structures can be used for the following:
  - Basic research – characterisation of the protein target (academic and industrial)
  - Lead discovery – identification of compounds (industrial and academic)

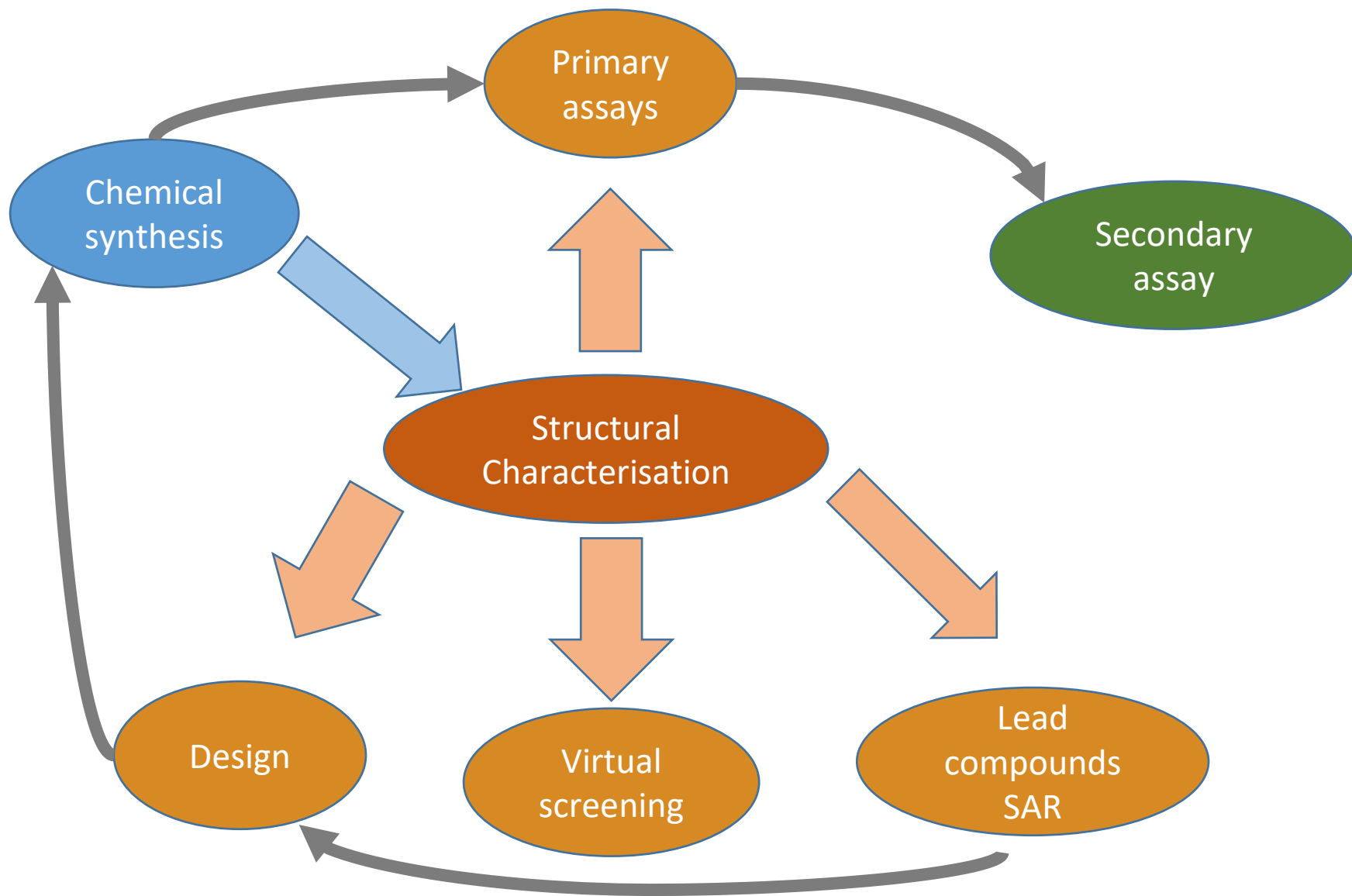

Structure is integral to many steps of the discovery process

# How structures contribute to screening



## Structural aided drug design

- Crystal structures help design molecules

## Virtual screen

- Docking models: use of a virtual compound library with the X-ray structure of the protein
- If have a known ligand, as a base to develop further compounds on

## Fragment screen



- Soak small compounds into crystals with often low mM activity

- Bound ligand in structure used to help predict where modifications could be added to provide increased potency or selectivity
- Often used as an adjunct to other screening strategies within big pharma
- Can provide the starting structures for a focused screen without the need to use expensive large library screens
- Can also be used to look for novel patent space around existing compound structures
- Identify new binding sites
- Identify novel chemistry for known sites
- Join selected fragments together to fit into the chemical space to increase potency


# SAR – Structure activity relationships




- Link between the drug target and the lead compound.
- Structures of the target with and without the lead compound act as a guide that can allow computational chemists and molecular modellers to identify binding site interactions

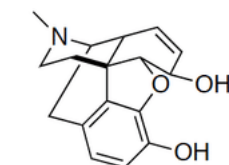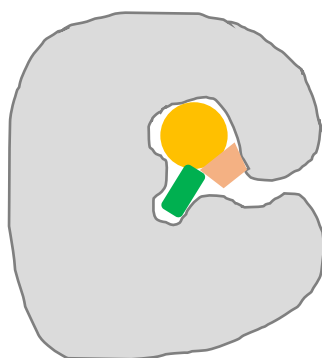


morphine  
relative potency = 1  
addictive

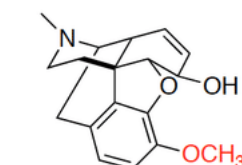


codeine  
relative potency = 0.2  
addictive

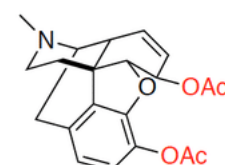


heroin  
relative potency = 2  
addictive

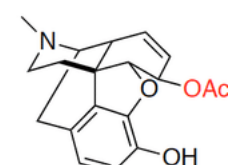


acetylmorphine  
relative potency = 4  
addictive

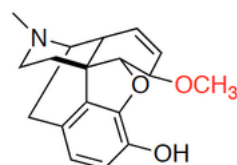


heterocodeine  
relative potency = 5  
addictive

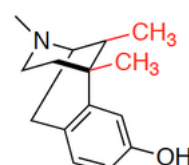

metazocine  
relative potency = 1  
addictive

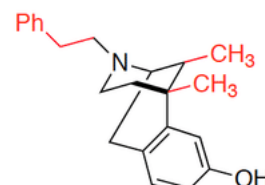

phenazocine  
relative potency = 4  
non-addictive

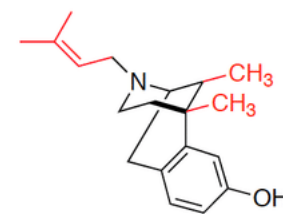


phenazocine  
relative potency = 0.33  
non-addictive




## What are Biologics?

A wide range of products such as vaccines, blood and blood components, allergenics, somatic cells, gene therapy, tissues, and recombinant therapeutic protein

While structure is not necessary for the development of biologics it is a useful analytical tool.

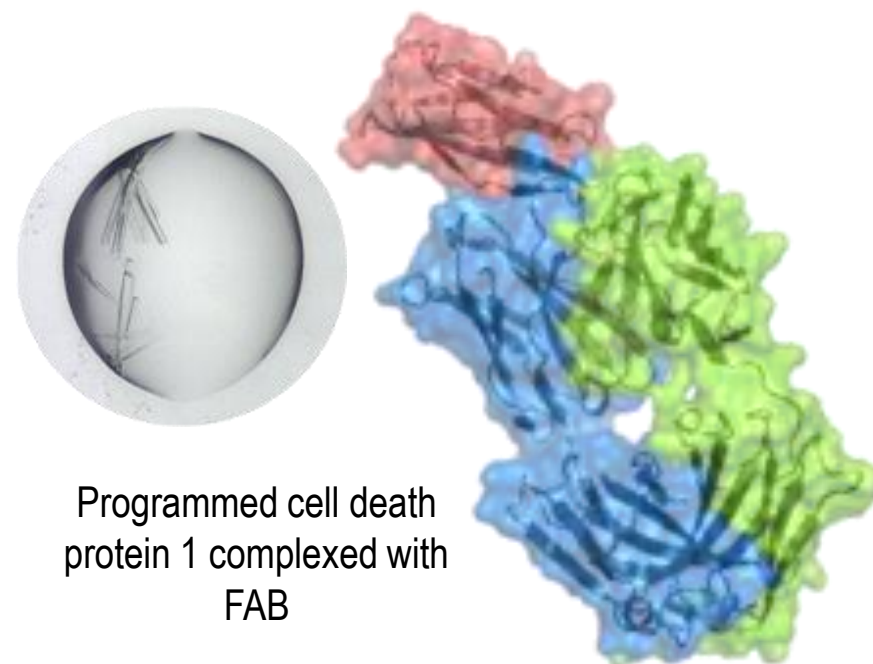

- Programmed Cell Death Protein 1 is a receptor protein on the surface of T cells
- Involved in the suppression of the immune system
- Blocking PD1 has been a successful approach to enhance anti-cancer immune responses
- Also interest in developing molecules that can trigger PD1 signalling to treat autoimmune diseases

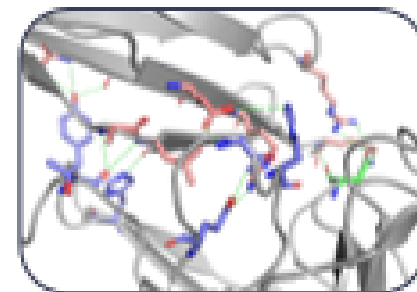



Crysalin Ltd

Venture capital funded spin out from University of Oxford



Founded by John Sinclair and Martin Noble

- Aim: to provide rapid determination of protein structures, irrespective of target class, at sub 3Å resolution


# How to achieve this?

Crysalins are a biological nanomaterial – a self assembling protein lattice

- Remove need for protein crystallisation
- No reproducibility issues (poor/variable resolution, etc)
- No target class restriction (ion channels/GPCRs)
- Remove need for large quantities of highly purified protein
- Increase speed of process
- Increase certainty of result


***An industrially robust and consistent process with predictable results***

***The freedom to choose which target to work on due to scientific validity of the target not the technical limitations of the methodology***

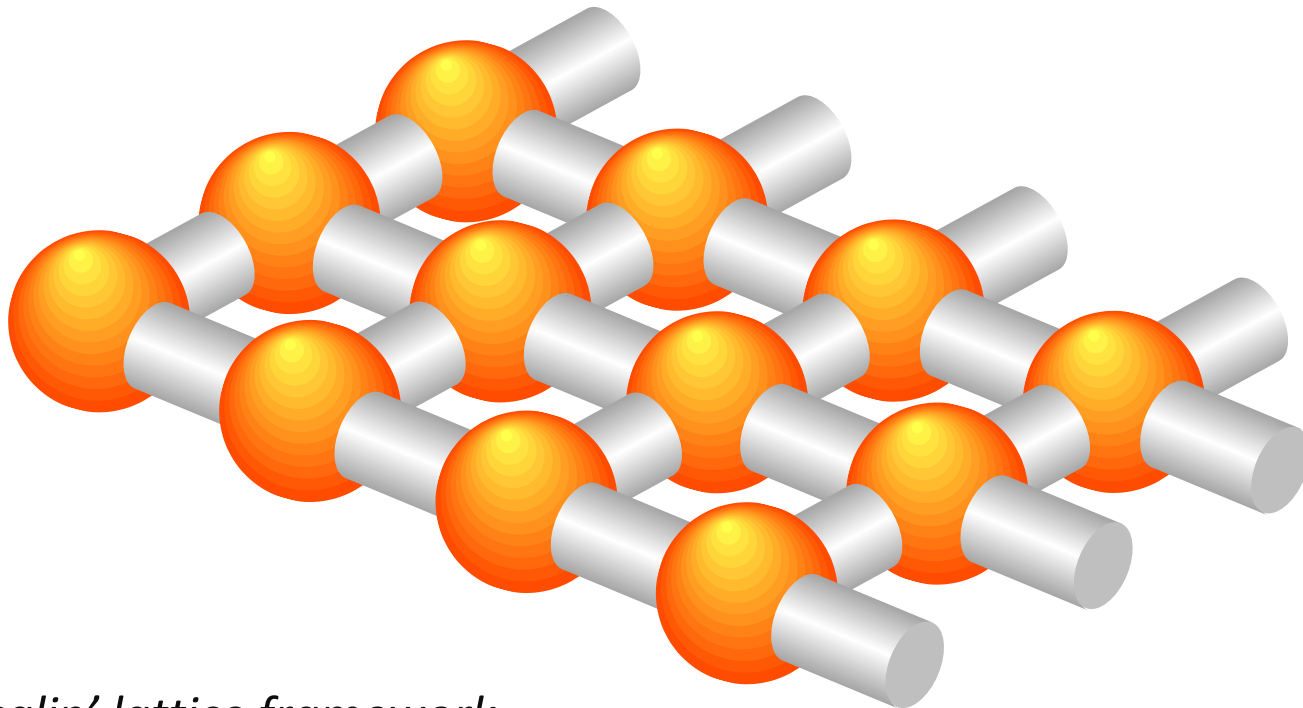


*The 'Crysalin' lattice framework...*

# Technical concept 2D

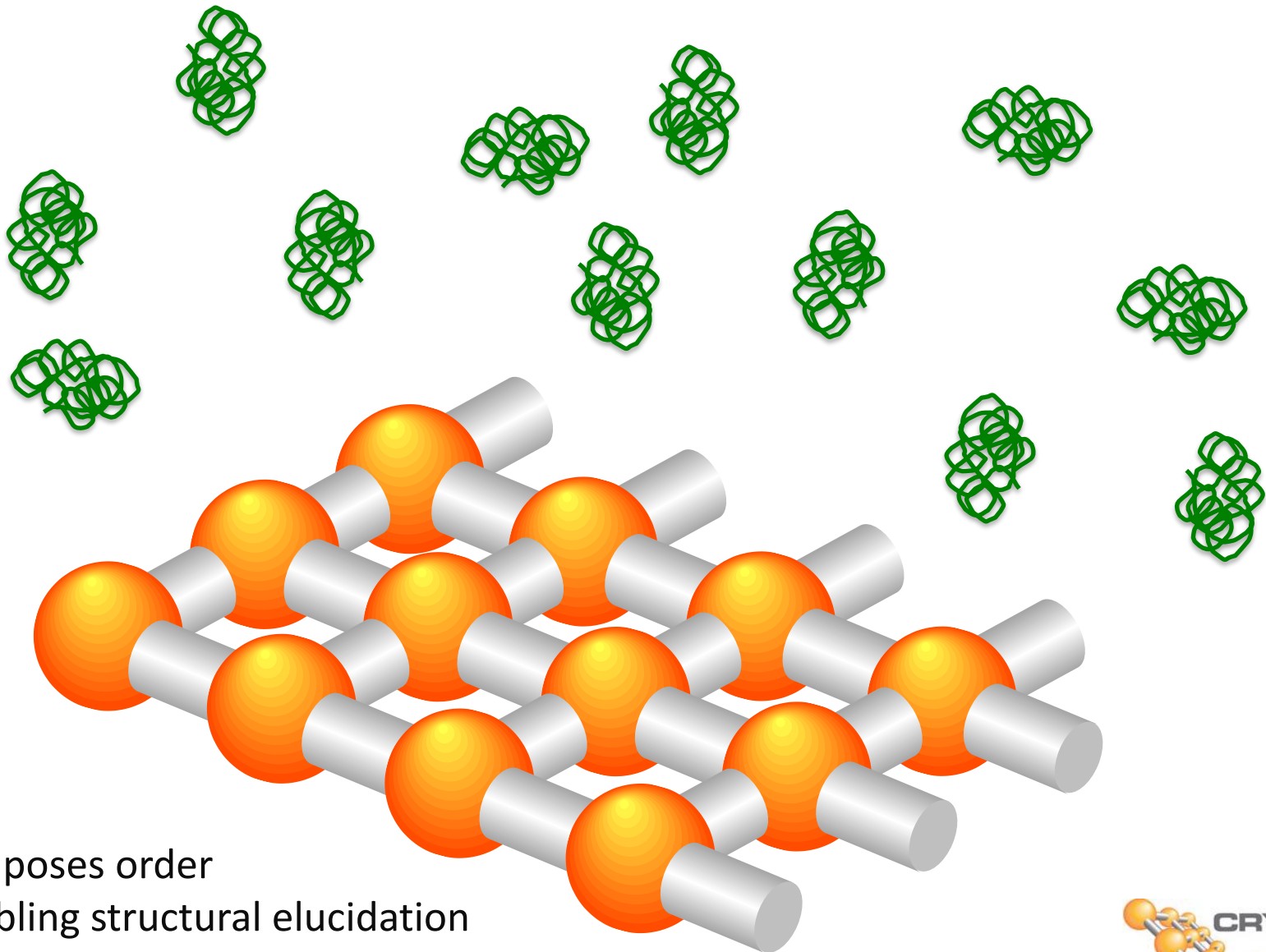


# 2D Crysalin lattice

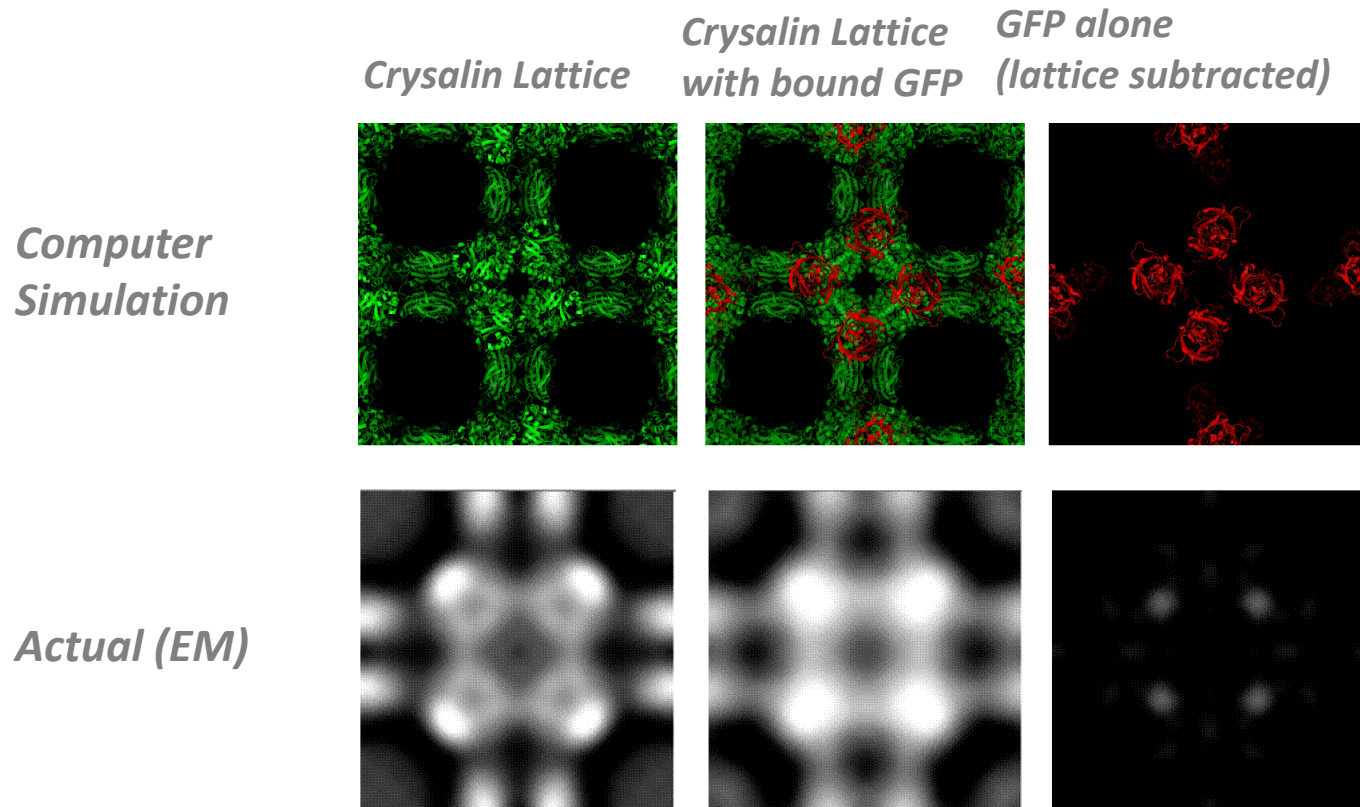

***Crysalin lattice manufacture sufficient to generate medium resolution structures (circa 12Å)\****

\* *Generation of protein lattices by fusing proteins with matching rotational symmetry*, John C. Sinclair, Karen M. Davies, Catherine Vénien-Bryan & Martin E. M. Noble, Nature Nanotechnology 6, 558–562 (2011)


# Technical concept 3D

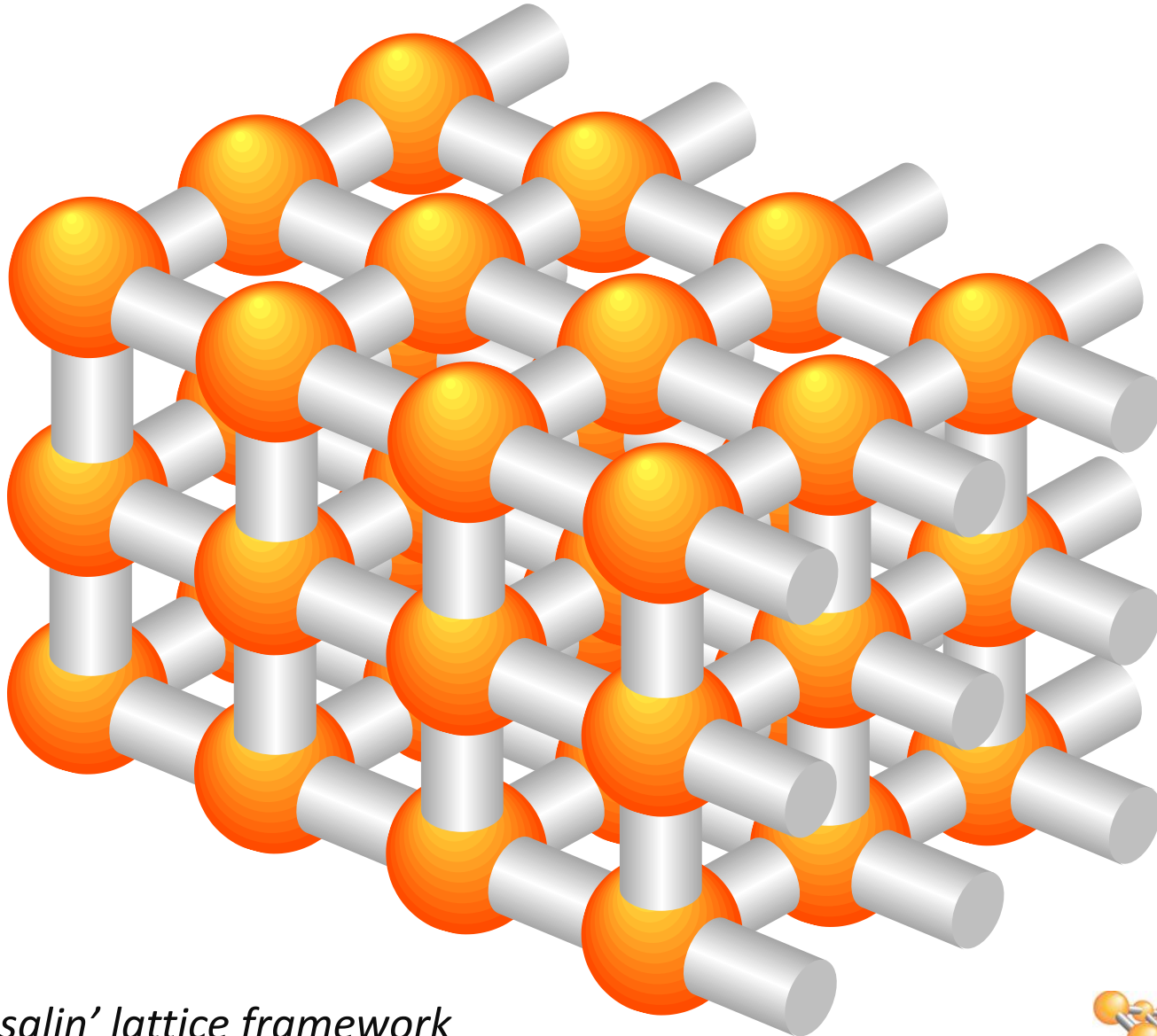




*The 3D 'Crysalin' lattice framework*




# Crysalin formation

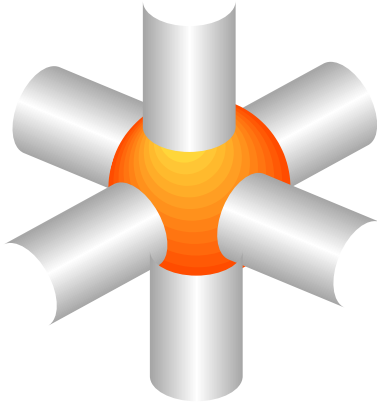

- Make the lattice components
- Crysalin lattice self assembles under the correct conditions
- Crysalins form as crystalline cubes varying in size from  $20\mu\text{m}$  to  $300\mu\text{m}$
- Typical diffraction ranges of base lattice from  $3.5\text{-}2.7\text{\AA}$  (best diffraction to  $2.5\text{\AA}$ )

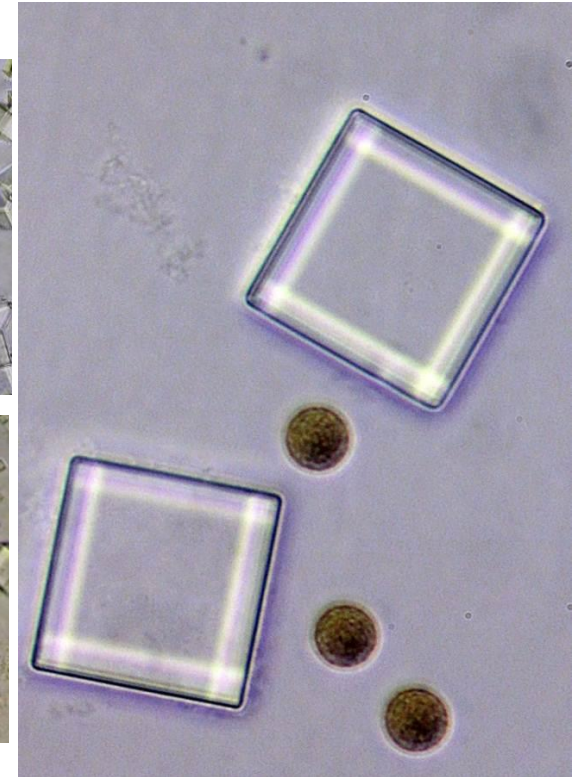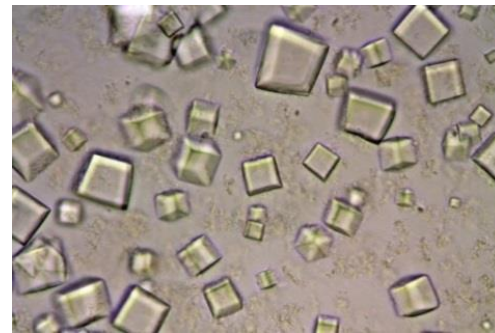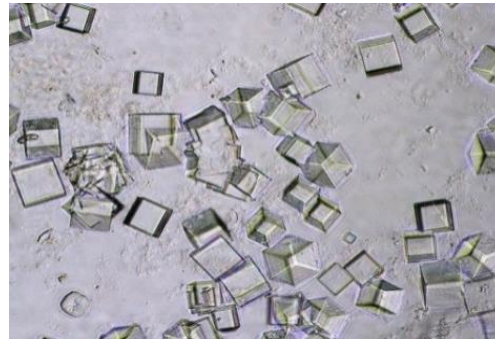

# XChem Platform

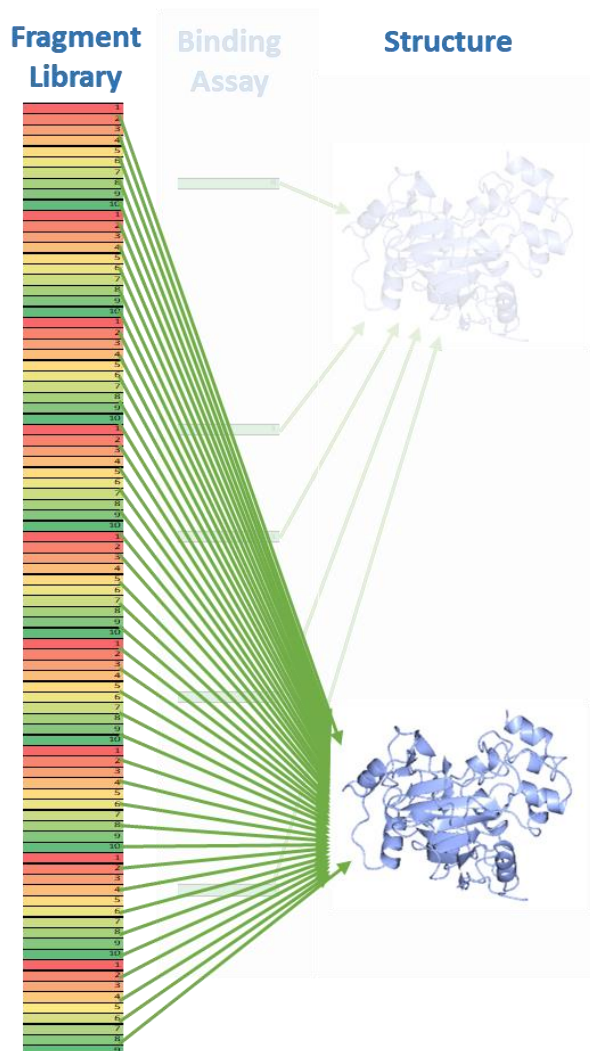

## Standard practice

- Cascade of biophysical methods
  - SPR, NMR, MST...
- Require crystal structures

## Via crystallography

- Very sensitive method
- High compound concentration
- Significant experimental overheads

## XChem platform


- >100 campaigns completed
  - 2-20% hit rate (project dependent)

# Overview of process


XChem Lab






Beamline I04-1

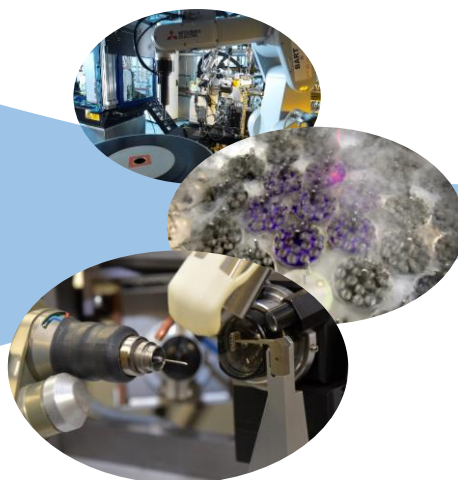

Example readout

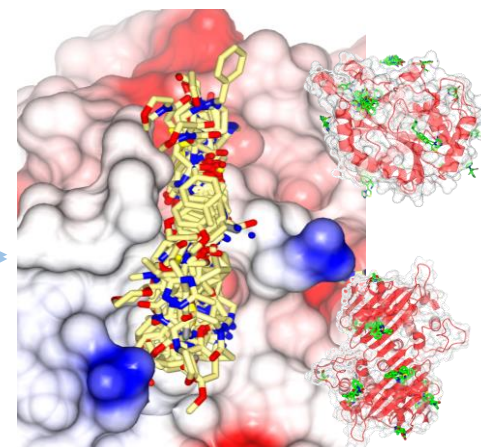

DSiP library in 1536 plate

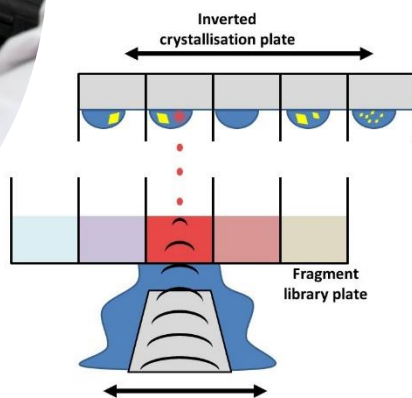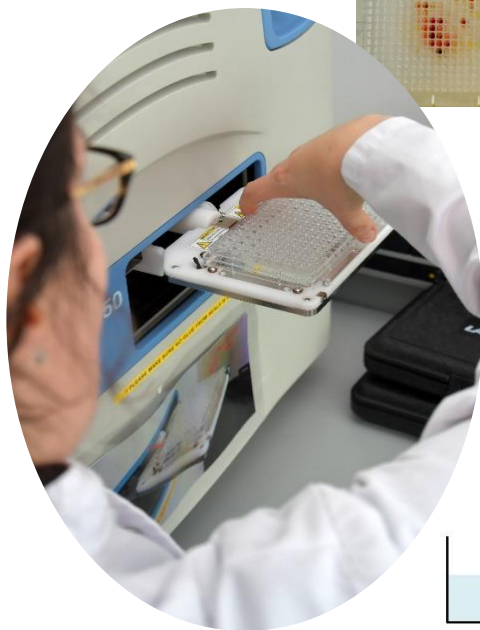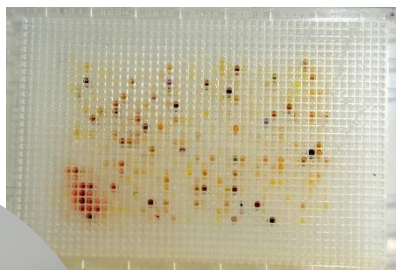


- Libraries
  - DSiP
    - ~950 compounds
    - 500mM in DMSO
  - External libraries
  - Cherry pick compounds
- Echo (Labcyte)
  - ~10min to transfer our DSiP library






- Crystal SHIFTER
  - Motorised X&Y stage under a microscope
  - Allows for rapid crystal harvesting
  - Touch screen graphical interface
    - Record events (time stamps)
      - Crystal mounted
      - Crystal status
      - Compound status
- Linked to SoakDB
  - XChem data management system

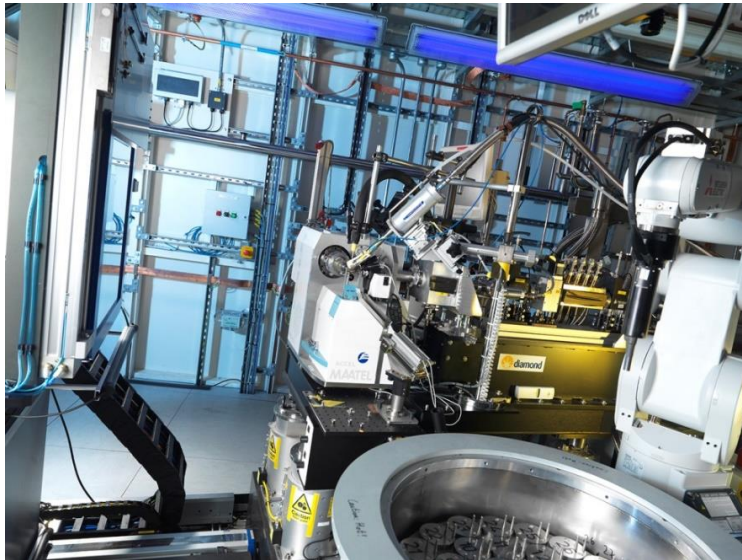


- I04-1 MX beamline
  - Automated loop centring
    - Optical or X-ray
  - Bart robot
    - Dewar capacity ~600 samples
  - Unattended data collection
  - Process up to 32 samples per hour
  - Eiger detector
    - 15s or 60s data collection
  - Auto-processing pipelines
    - Xia2, Dials and AutoProc



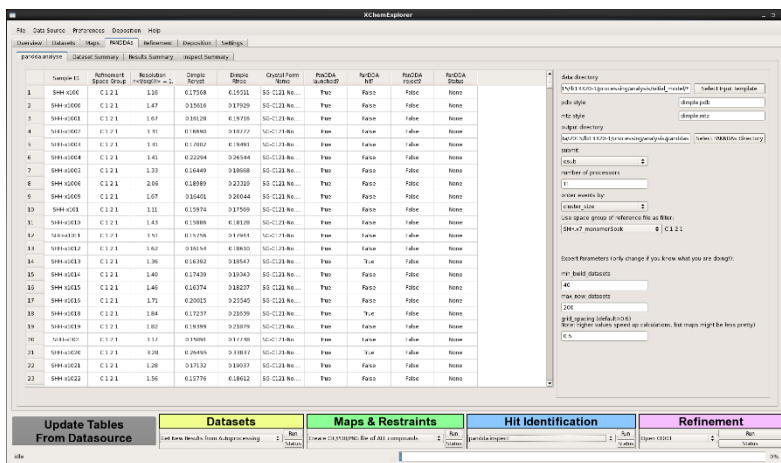

Sample ID	Substrate	Resolution (Å)	Crystals	Crystals (µm)	Crystals (nm)	PanDDA	PanDDA	PanDDA	PanDDA
						Hit	Hit	Hit	Hit
1	SHH-0001	1.15	1.15	1.15	1.15	True	False	False	False
2	SHH-0002	1.15	1.15	1.15	1.15	True	False	False	False
3	SHH-0003	1.15	1.15	1.15	1.15	True	False	False	False
4	SHH-0004	1.15	1.15	1.15	1.15	True	False	False	False
5	SHH-0005	1.15	1.15	1.15	1.15	True	False	False	False
6	SHH-0006	1.15	1.15	1.15	1.15	True	False	False	False
7	SHH-0007	1.15	1.15	1.15	1.15	True	False	False	False
8	SHH-0008	1.15	1.15	1.15	1.15	True	False	False	False
9	SHH-0009	1.15	1.15	1.15	1.15	True	False	False	False
10	SHH-0010	1.15	1.15	1.15	1.15	True	False	False	False
11	SHH-0011	1.15	1.15	1.15	1.15	True	False	False	False
12	SHH-0012	1.15	1.15	1.15	1.15	True	False	False	False
13	SHH-0013	1.15	1.15	1.15	1.15	True	False	False	False
14	SHH-0014	1.15	1.15	1.15	1.15	True	False	False	False
15	SHH-0015	1.15	1.15	1.15	1.15	True	False	False	False
16	SHH-0016	1.15	1.15	1.15	1.15	True	False	False	False
17	SHH-0017	1.15	1.15	1.15	1.15	True	False	False	False
18	SHH-0018	1.15	1.15	1.15	1.15	True	False	False	False
19	SHH-0019	1.15	1.15	1.15	1.15	True	False	False	False
20	SHH-0020	1.15	1.15	1.15	1.15	True	False	False	False
21	SHH-0021	1.15	1.15	1.15	1.15	True	False	False	False
22	SHH-0022	1.15	1.15	1.15	1.15	True	False	False	False

Krojer *et al*, 2017, ActaD

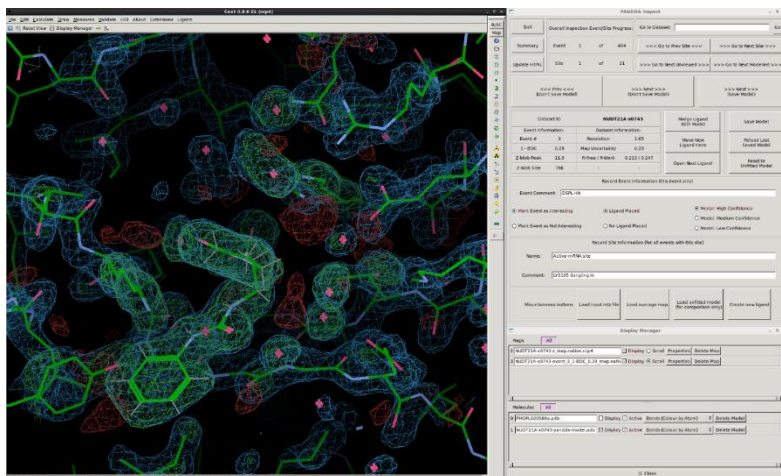

PanDDA. Pearce *et al*, 2017, Nat. Comm.

## • XChemExplorer

### • Record all information

- Same SQLite database as SoakDB

### • Rapid review of data collection

- Selection of processed files

### • Generation of maps

- Run Dimple on user's PDB model

### • PanDDA user interface


#### • Pre-PanDDA

- Generate accurate PDB model of "background state"


#### • Coot Plugin

- Facilitate ligand fitting
- Quickly browse through events and record comments

## Results from JMJD2 campaign (SGC)

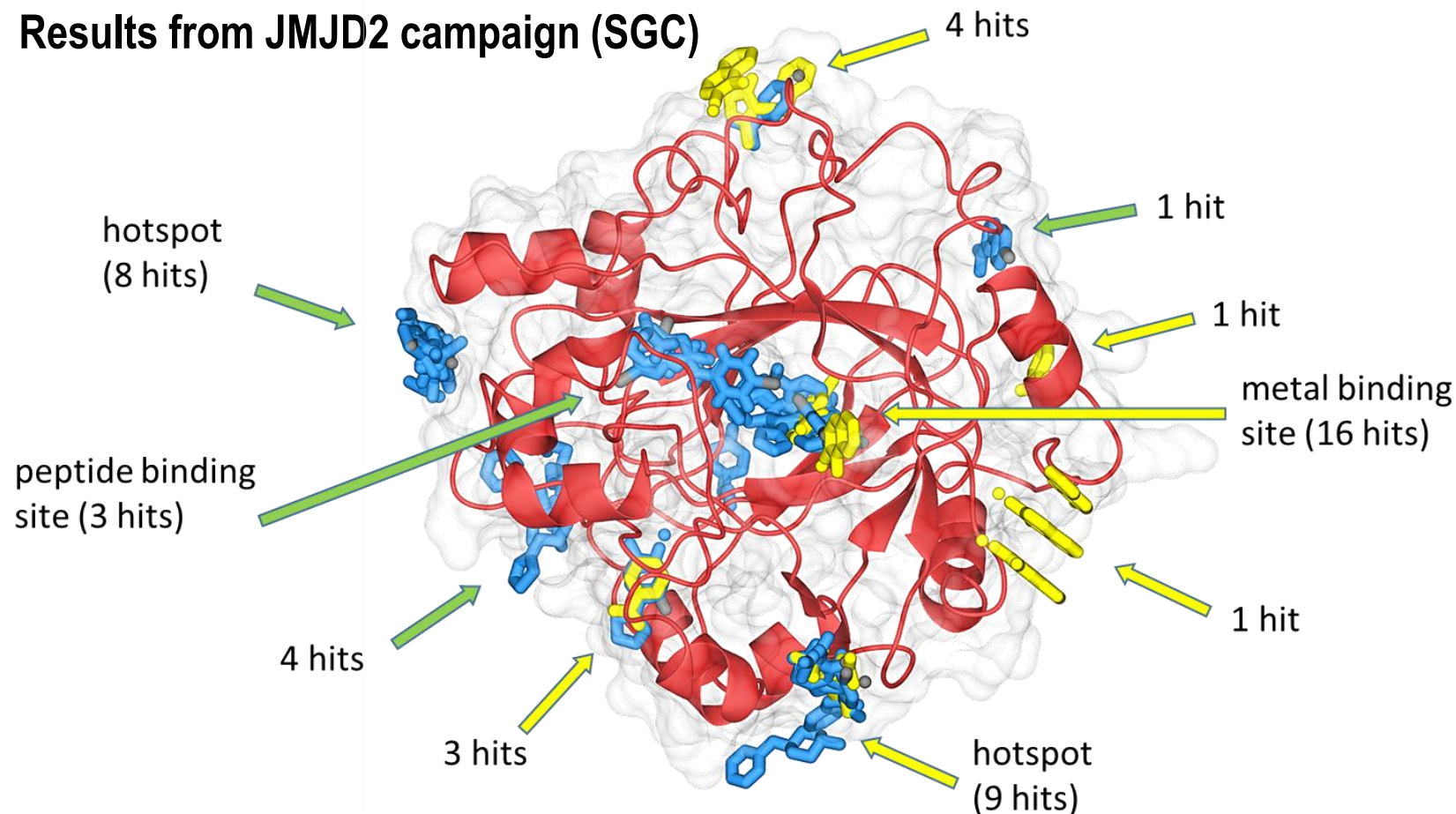

# Order of magnitude speedup

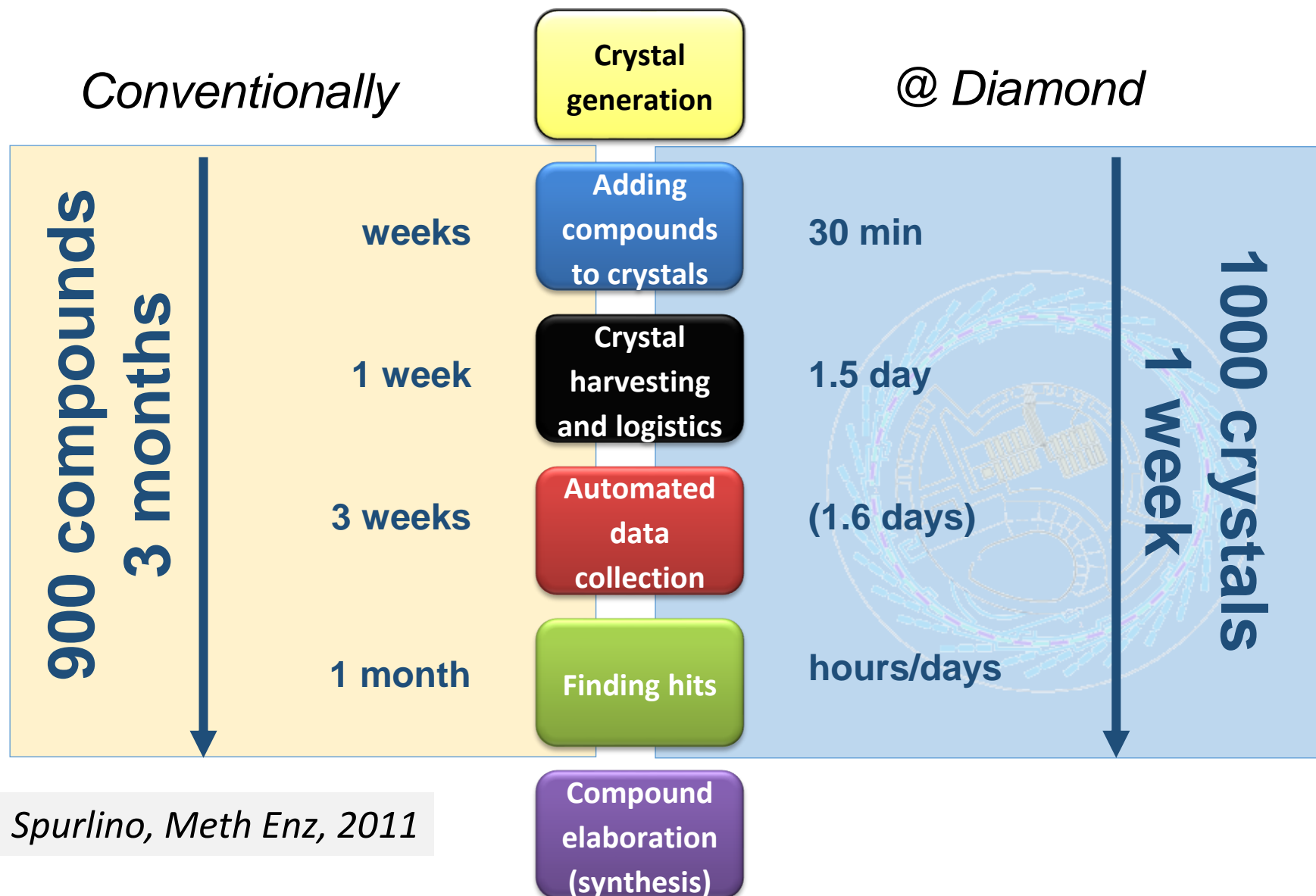

# How is the platform used by industry?



- Trained on the platform and use it with support in place
- Send projects for “full service” access
  - Project run by members of Diamond’s ILO team
  - This often includes data analysis
- Companies run multiple campaigns per year
- Incorporating into their general screening protocols
- Frequently have follow up compounds sent for testing
- Over subscribed – often running 12-14 campaigns at different stages at any one time
- Increasing enquires from virtual and AI companies

- Rapid turnover of data
- Different approach to data collection
- Quick analysis of data to report to client/manager/board
- Steer the scientific direction
- Can be working on a well defined system, but also developing that system
- Expectation of success
- If a project isn't working, will be stopped
- Interactions across teams
- Rapid problem solving
- Working on multiple projects
- Development of transferable skills, eg.
  - Management
  - Business development

Questions?