

Refinement

Garib N Murshudov

MRC-LMB

Cambridge


Contents

- 1) Purpose of refinement**
- 2) Model and Data**
- 3) Model: what do we know about macromolecules**
- 4) Bulk solvent modelling**
- 5) Crystallographic data and their property**
- 6) Model parameters**

Purpose of refinement


Crystallographic refinement has two purposes:

- 1) Fit chemically and structurally sensible atomic model into observed– X-ray crystallographic data**
- 2) To calculate best possible (electron density) map so that atom model can be criticised**

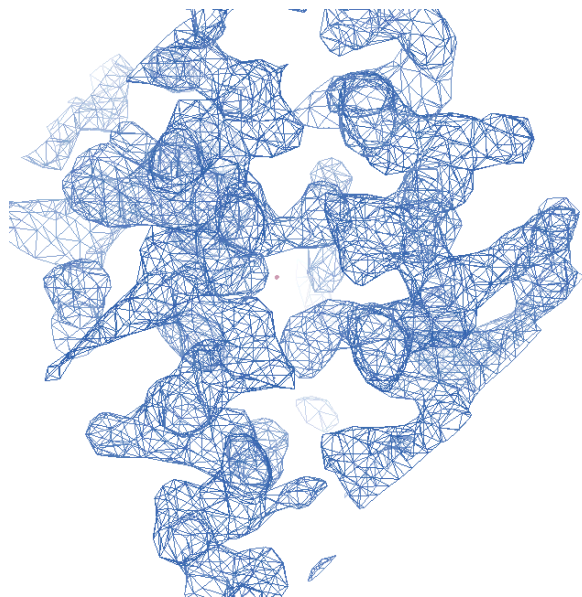


Data

+

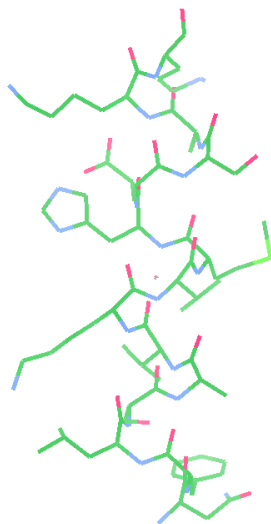

Atomic model

=

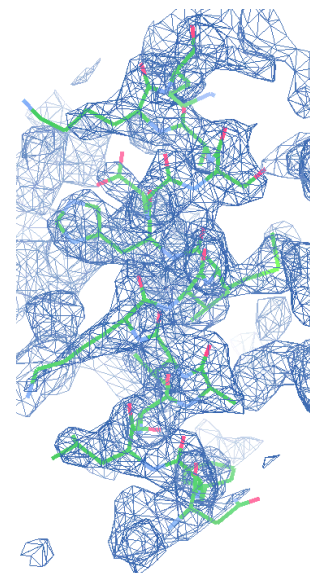

Fit and refine

We want to fit currently available model to the data and calculate differences between them.


To do this fit properly we must use as much as possible information about model and data.

What do we know about macromolecules?

- 1) Macromolecular consists of atoms that are bonded to each other in a specific way**
- 2) If there are two molecules with sufficiently high sequence identity then it is likely that they will be similar to each other in 3D**
- 3) It is highly likely that if there are two copies of a the same molecule they will be similar to each other (at least locally)**
- 4) Oscillation of atoms close to each other in 3D cannot be dramatically different**
- 5) Proteins tend to form secondary structures**
- 6) DNA/RNA tend to form base-pairs, stacked bases tend to be parallel**

Why restraints: Chemical knowledge

- Distance between atoms 1.3Å. B values 20 and 50

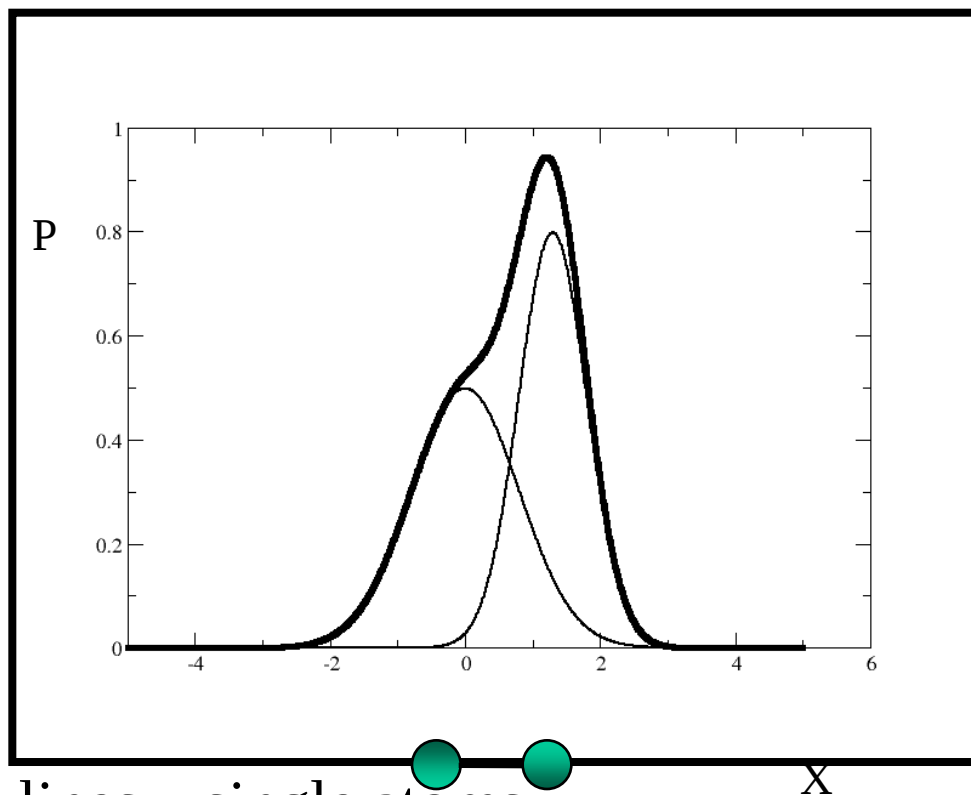

- Thin lines – single atoms
- Bold line - sum of the two atoms

NCS

Three ways of dealing with NCS

- 1) NCS constraints: copies of molecules are considered to be exactly same. Only one set of atomic parameters per molecule is refined, other copies are kept to be exactly same
- 2) NCS restraints: Molecules are superimposed and difference between corresponding atoms after superposition minimised.
- 3) NCS local restraints: Molecules are assumed to be locally similar, globally they may be different

Auto NCS: local and global


1. Align all chains with all chains using Needleman-Wunsh method
2. If alignment score is higher than predefined (e.g.80%) value then consider them as similar
- 3.Find local RMS and if average local RMS is less than predefined value then consider them aligned
4. Find correspondence between atoms
5. If global restraints (i.e. restraints based on RMS between atoms of aligned chains) then identify domains
- 6.For local NCS make the list of corresponding interatomic distances (remove bond and angle related atom pairs)
- 7.Design weights

The list of interatomic distance pairs is calculated at every cycle

Auto NCS

Global RMS is calculated using all aligned atoms.

Local RMS is calculated using k (default is 5) residue sliding windows and then averaging of the results

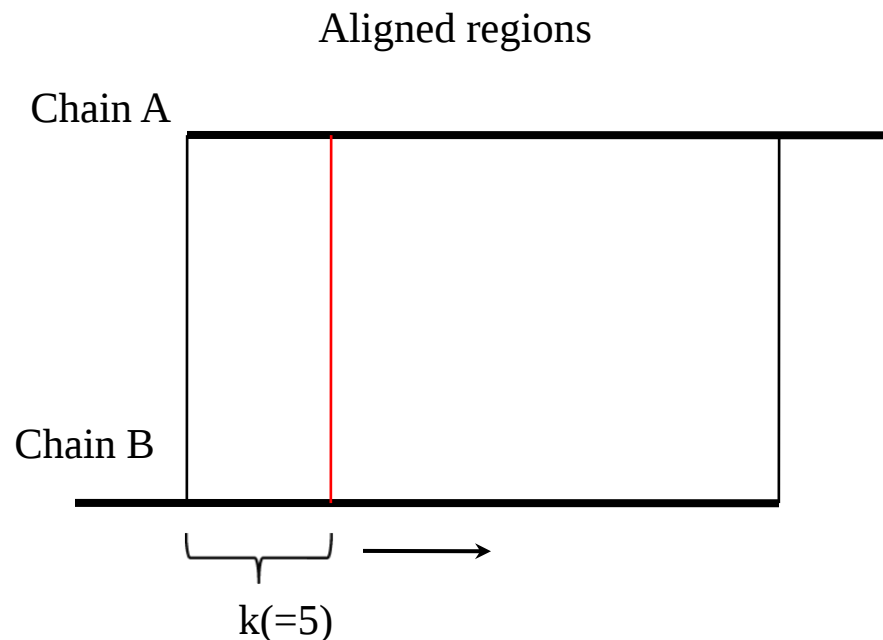

$$Ave(RmsLoc)_k = \frac{1}{N - k + 1} \sum_{i=1}^{N - k + 1} RmsLoc_i$$

$$RMS = Ave(RmsLoc)_N$$

Auto NCS: Iterative alignment

Example of alignment: 2vtu.

There are two chains similar to each other. There appears to be gene duplication

RMS – all aligned atoms

Ave(RmsLoc) – local RMS



***** Alignment results *****

: N:	Chain 1 :	Chain 2 :	No of aligned	:Score :	RMS	:Ave(RmsLoc):

: 1 :	J(131 - 256) :	J(3 - 128) :	126 :	1.0000 :	5.2409 :	1.6608 :
: 2 :	J(1 - 257) :	L(1 - 257) :	257 :	1.0000 :	4.8200 :	1.6694 :
: 3 :	J(131 - 256) :	L(3 - 128) :	126 :	1.0000 :	5.2092 :	1.6820 :
: 4 :	J(3 - 128) :	L(131 - 256) :	126 :	1.0000 :	3.0316 :	1.5414 :
: 5 :	L(131 - 256) :	L(3 - 128) :	126 :	1.0000 :	0.4515 :	0.0464 :

Auto NCS: Conformational changes

In many cases it could be expected that two or more copies of the same molecule will have (slightly) different conformation. For example if there is a domain movement then internal structures of domains will be same but between domains distances will be different in two copies of a molecule

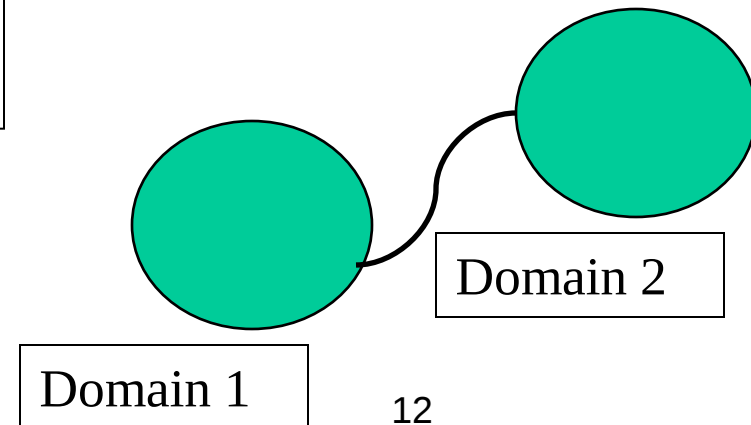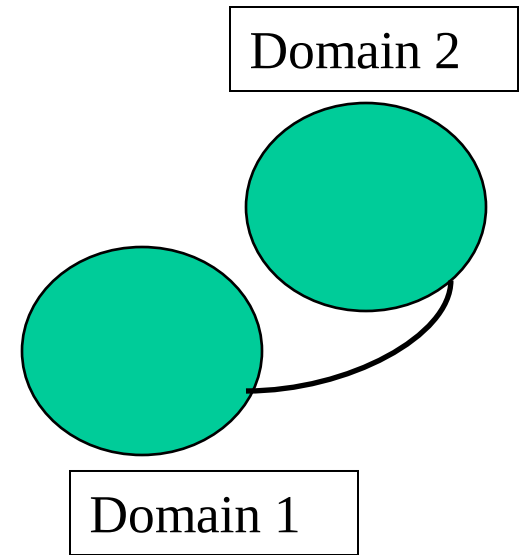

External (reference) structure restraints

Restraints to external structures are generated by the program ProSmart:

- 1) Aligns structure in the presence of conformational changes. Sequence is not used
- 2) Generates restraints for aligned atoms
- 3) Identifies secondary structures (at the moment helix and strand, but the approach is general and can be extended to any motif).
- 4) Generates restraints for secondary structures


Note 1: ProSmart has been written by Rob Nicholls and available from him and CCP4.

Note 2: Robust estimator functions are used for restraints. I.e. if differences between target and model is very large then their contributions are down-weighted

External Restraints

Restraint Parameter Estimation:


Atoms in homologous structure:

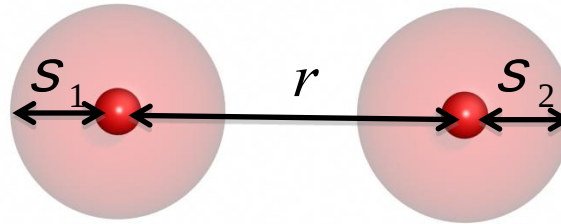

σ_i estimated as by Murshudov & Dodson (CCP4 Newsletter 1997)


Dependencies:

- B-value
- Resolution
- Model/data completeness
- Data quality

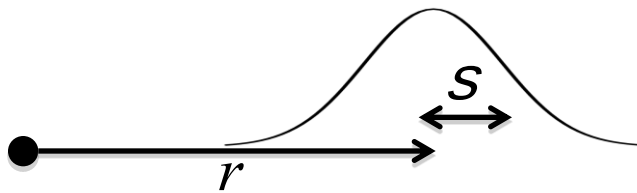

Plus a modification that increases robustness to outliers (use of Geman-McClure function)

External Restraints

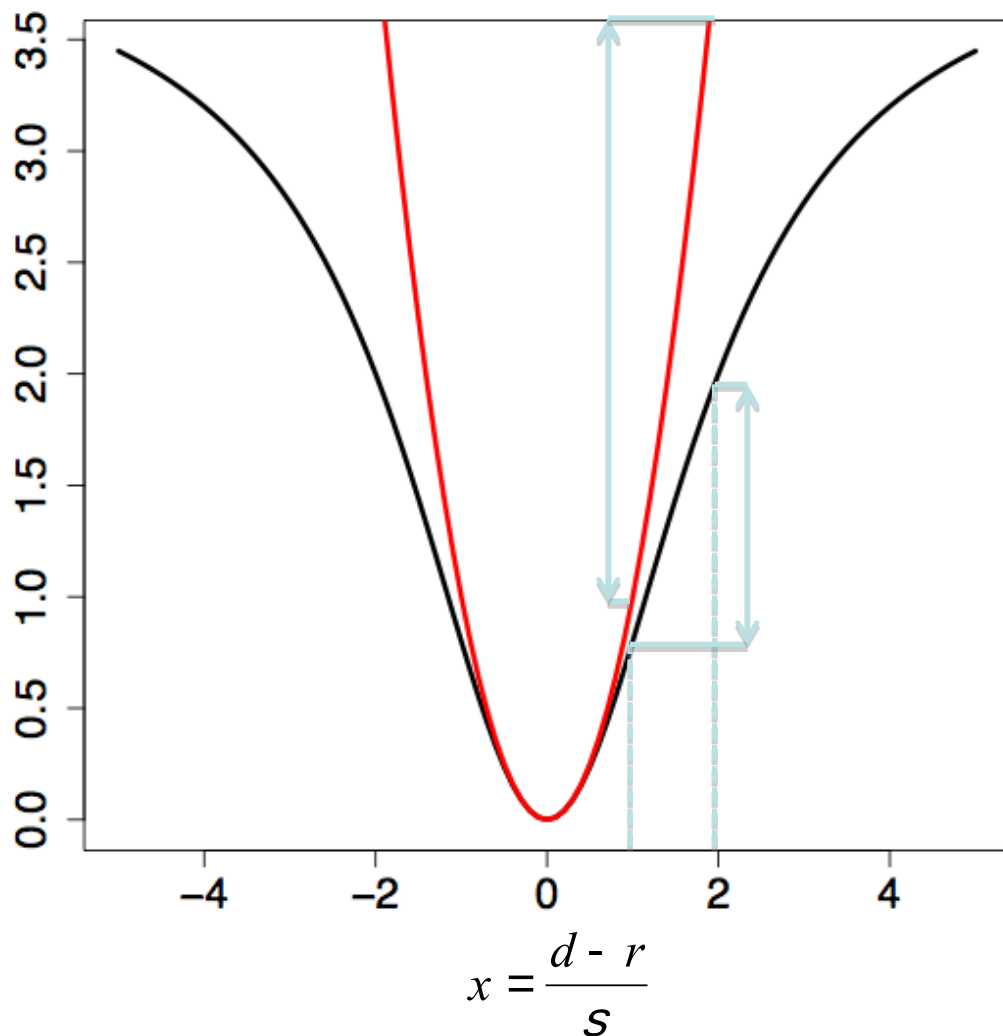


Least squares: x^2

Geman-McClure: $\frac{x^2}{1 + wx^2}$


Also use external B-value restraints and Van-der-Waals repulsions in recent versions of REFMAC5.

External Restraints

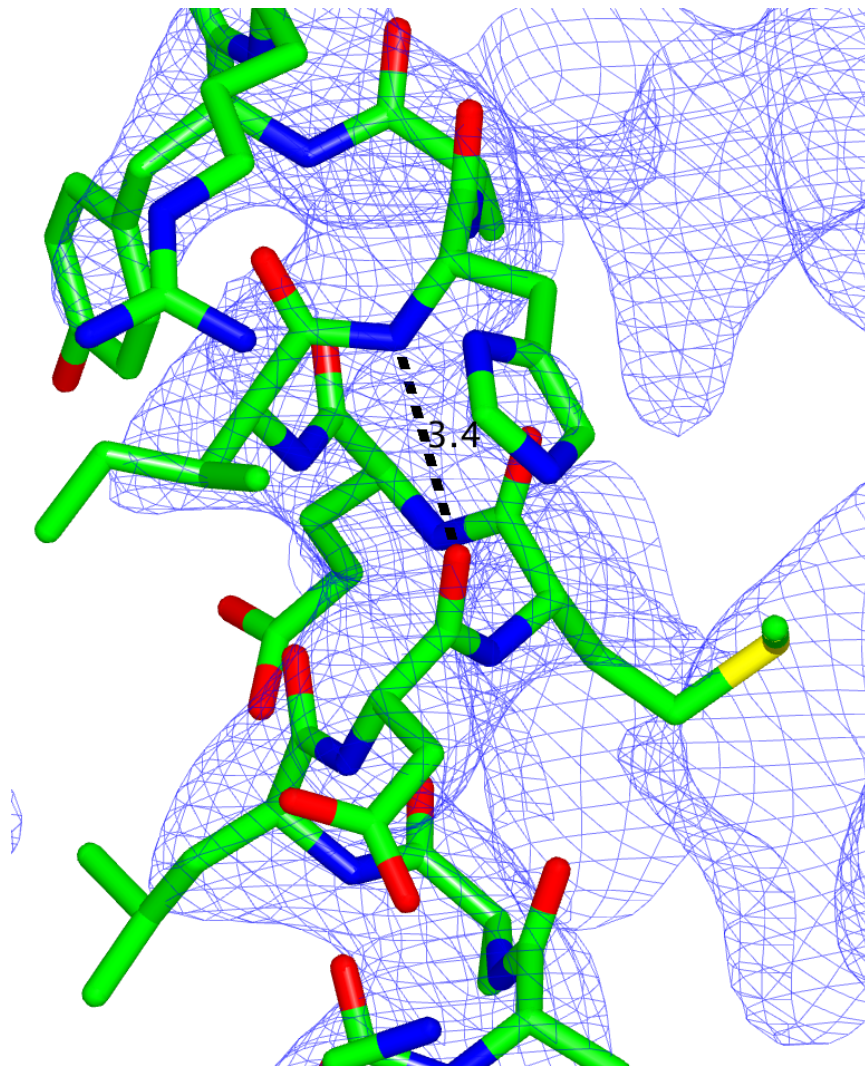

3g4w – 3.7 Å

Prior information:

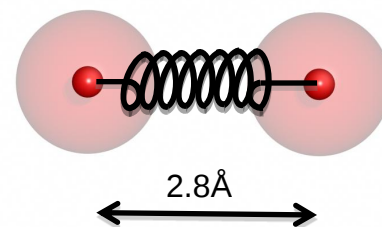


External Restraints

Example: Ovotransferrin



External Restraints

Example: Ovotransferrin

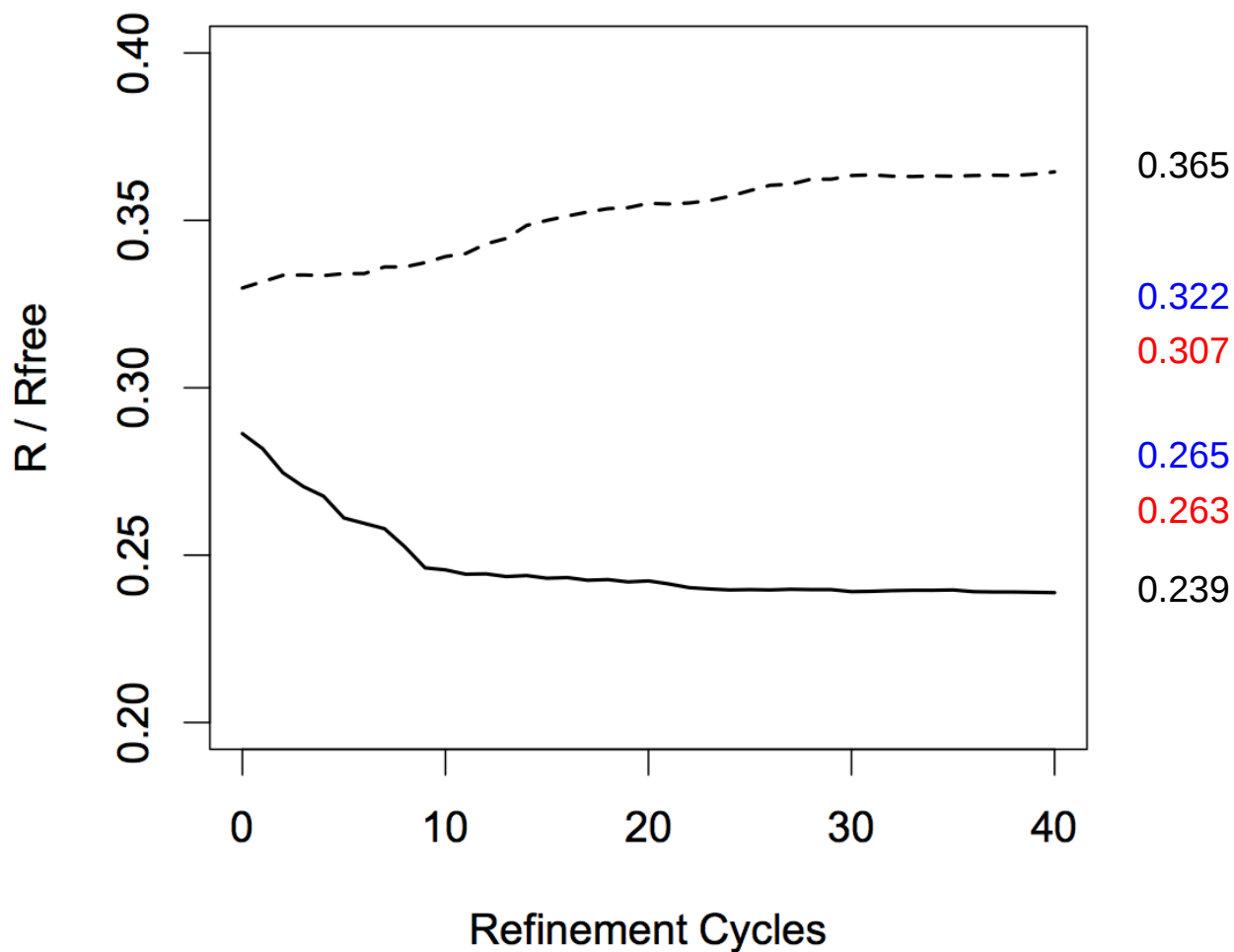


External Restraints

Example: Ovotransferrin


Original Structure

$R/R_{\text{free}} : 0.286/0.330$

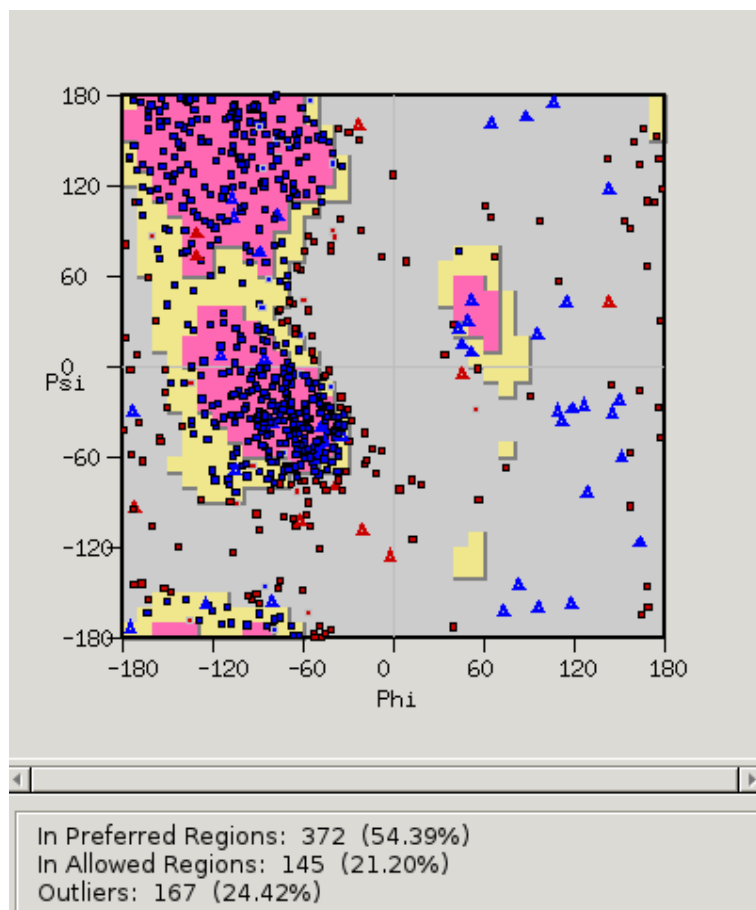

Re-refined with External Restraints


$R/R_{\text{free}} : 0.263/0.307$

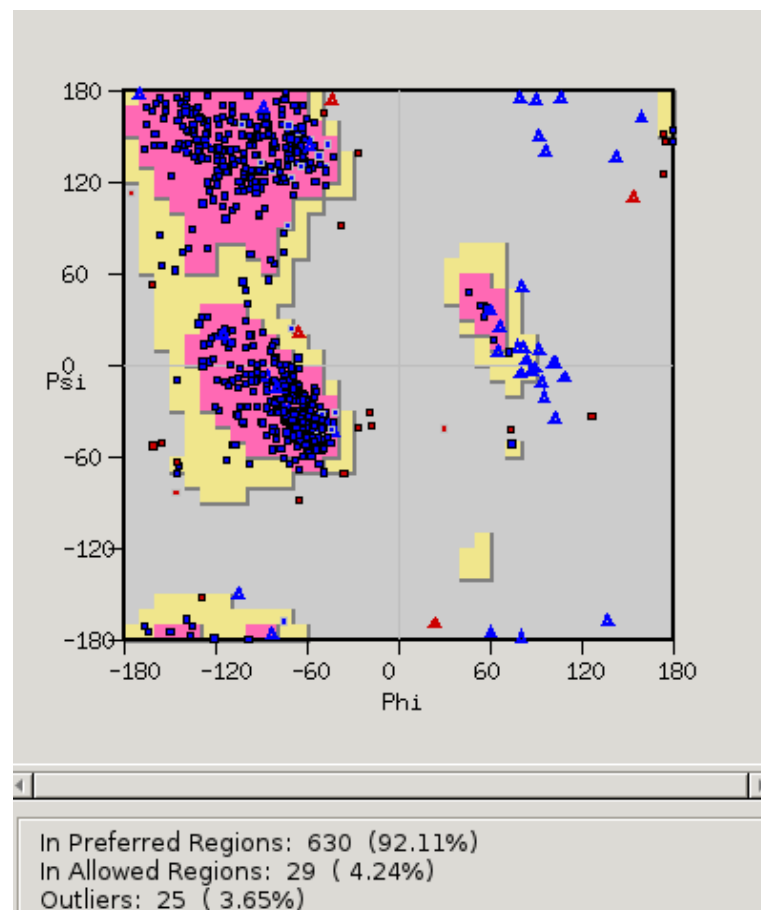

External Restraints


Original Structure

R/R_{free} : 0.286/0.330

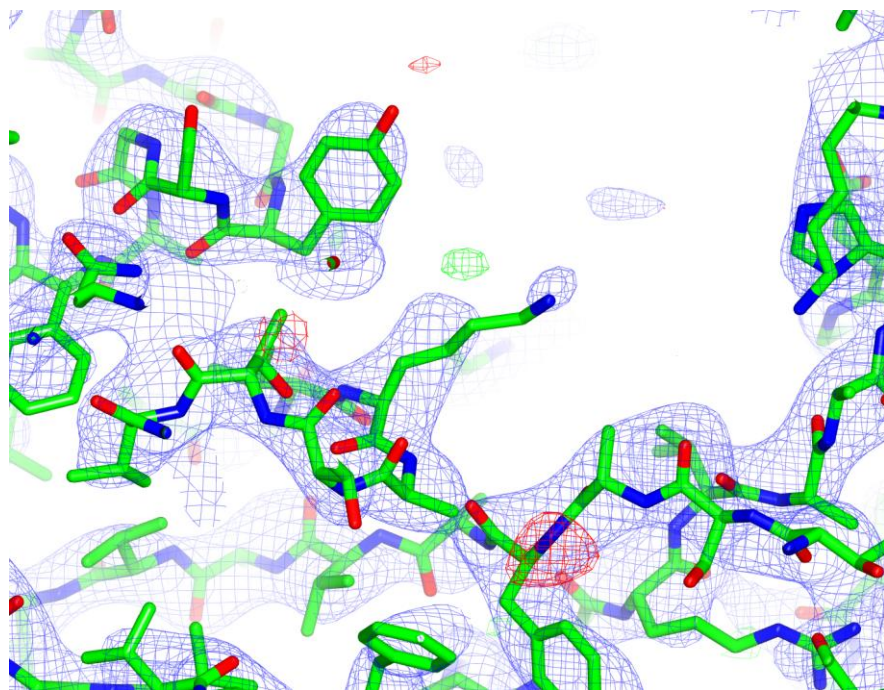


Re-refined Structure

R/R_{free} : 0.252/0.307

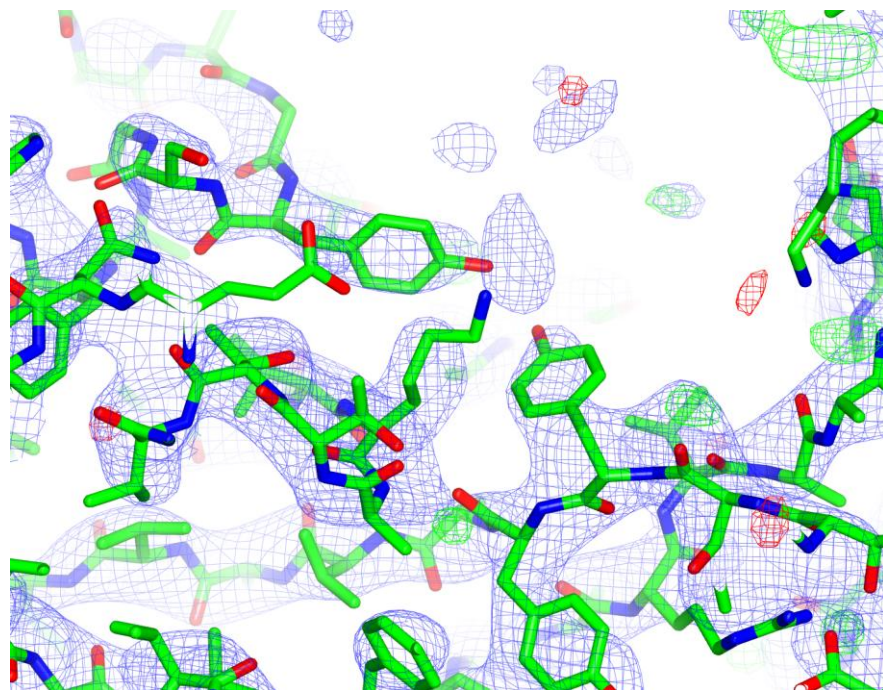

ProSMART Restraint Visualisation in Coot

Backbone
Restrains

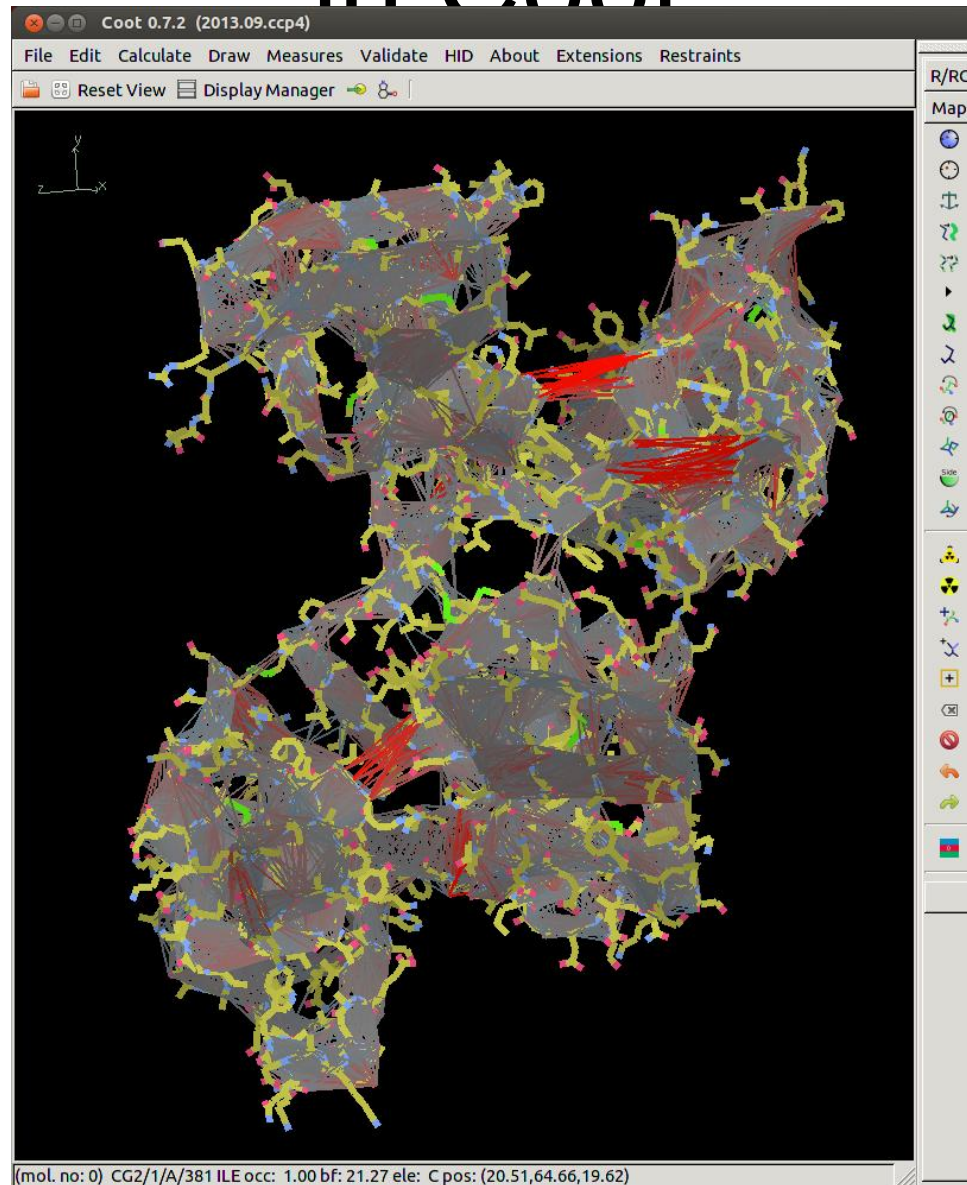

External Restraints

Summary:

Regularisation helps to stabilise refinement by reducing effective number of parameters, thus increasing the signal-to-noise ratio.

New techniques may:

- Reveal new map features
- Lead to model improvement

When refining at low resolution, check:

- Refinement statistics - Not always conclusive
- Geometry - Not always conclusive
- Electron density - Not always reliable

Conclusion: take care.

Basepair restraints

monomers C G label C:G

bond atom N4 C atom O6 G value 2.91 sigma 0.15

bond atom N3 C atom N1 G value 2.95 sigma 0.1

bond atom O2 C atom N2 G value 2.86 sigma 0.15


chiral atom N1 G atom C6 G atom C2 G atom N3 C value 0.0 sigma 0.5

chiral atom N3 C atom C4 C atom C2 C atom N1 G value 0.0 sigma 0.5

torsion atom C4 C atom N3 C atom N1 G atom C2 G value 180.0 sigma 10.0

Tools for generation of extra DNA/RNA restraints

- ❖ Generate extra restraints on sub-structures such as base-pairs, backbone conformation, stacking planes in DNA, RNA and their protein complex for refinement

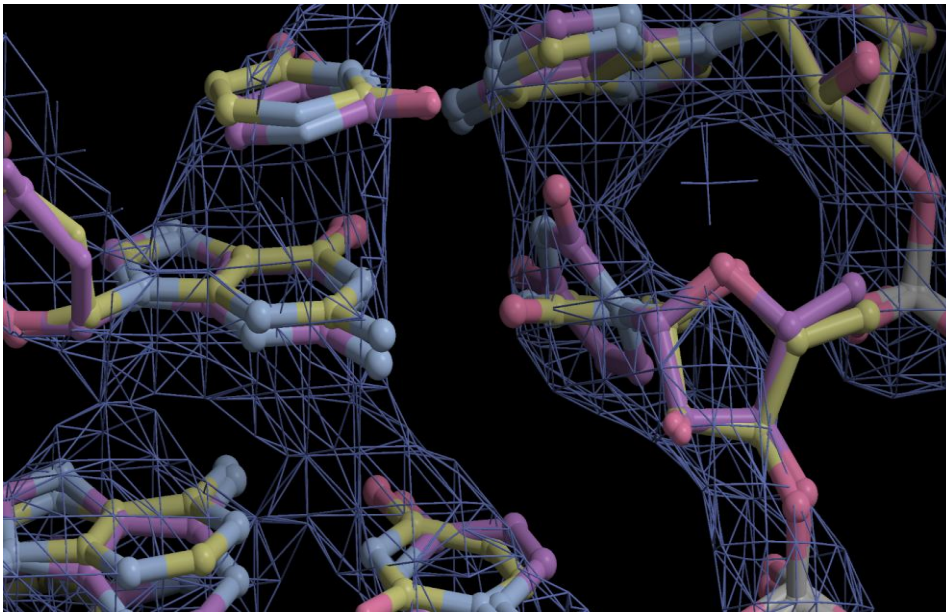


- ❖ Data from single-particle electron cryo-microscopy
- ❖ Resolution 3.2Å
- ❖ G:U base pair. Purple: before refinement. Yellow: after refinement. The base-pair is put to the same plane via refinement

Canonical and non-canonical base pairs can be restrained. In case of non-canonical base-pairing there might be several possible pairings. During refinement the closest pairing is selected.

Tools for generation of extra DNA/RNA restraints

- ❖ Generate extra restraints on sub-structures such as base-pairs, backbone conformation, stacking planes in DNA, RNA and their protein complex for refinement

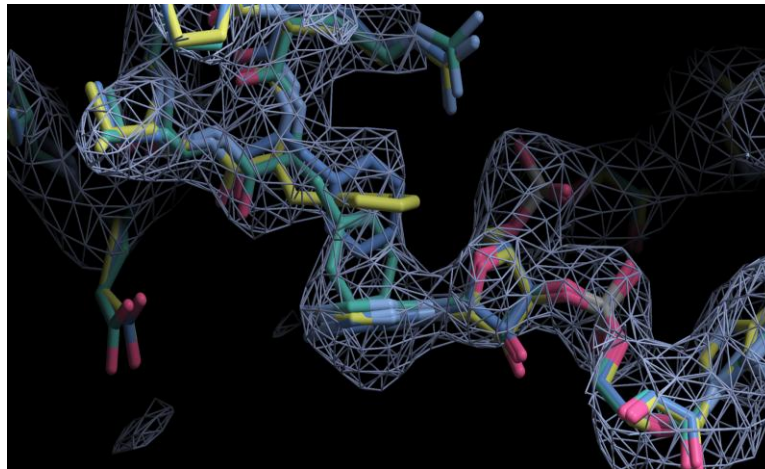

- ❖ Data from single-particle electron cryo-microscopy
- ❖ Resolution 3.2Å
- ❖ Green: before refinement. Blue: after refinement without stacking. Yellow: after refinement with stacking

Planes can be restrained to be parallel or to form a particular angle

Restraints to current distances (jelly-body)

The term is added to the target function:


$$\sum_{pairs} w(|d| - |d_{current}|)^2$$

Summation is over all pairs in the same chain and within given distance (default 4.2Å). $d_{current}$ is recalculated at every cycle. This function does not contribute to gradients. It only contributes to the second derivative matrix.


It is equivalent to adding springs between atom pairs. During refinement inter-atomic distances are not changed very much. If all pairs would be used and weights would be very large then it would be equivalent to rigid body refinement.

It could be called “implicit normal modes”, “soft” body or “jelly” body refinement.

Effect of “jelly” body refinement: Example is provided by A. Lebedev

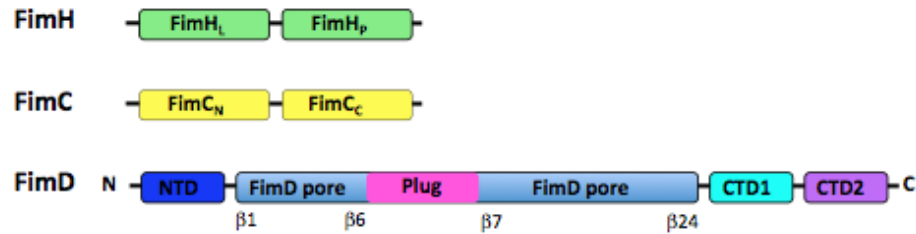



- Asymmetric unit two
- Resolution 2.8 Å

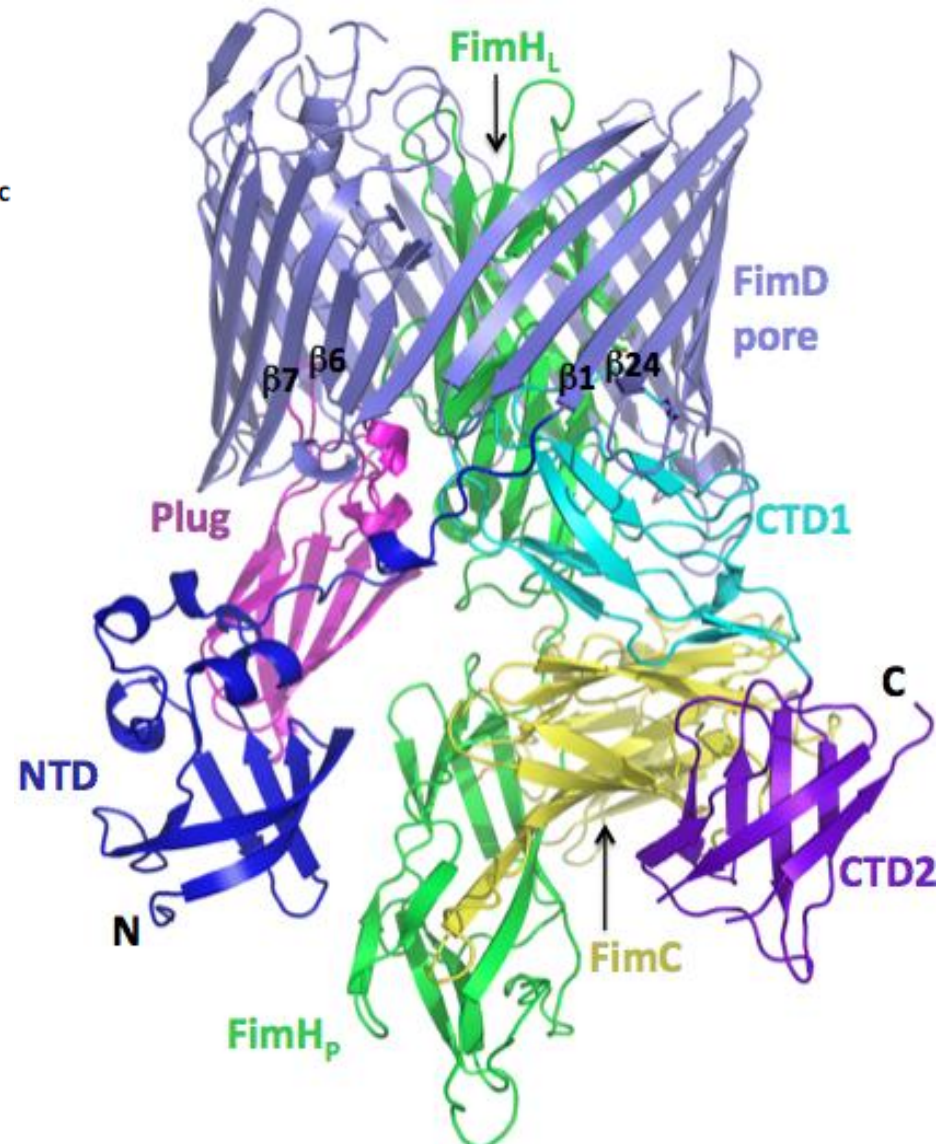

Phane et. al (2011) Nature, 474, 50-53

Usher complex structure solution

- 1. Conventional MR
 - FimC-N + FimC-C
 - FimH-L + FimH-P
 - FimD-Pore
- 2. Jelly body refinement (Refmac)
 - FimD-Pore

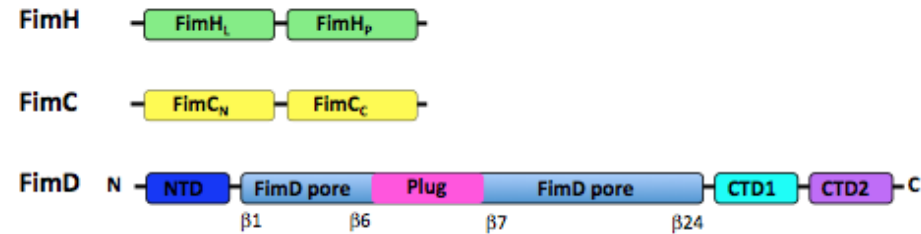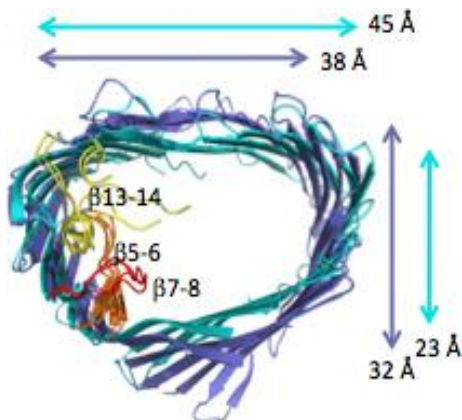

- 3. Fitting into the electron density
 - FimD-Plug
 - FimD-NTD
 - FimD-CTD-2
- 4. Manual building
 - FimD-CTD-1

Data and their property

- 1) **Amplitudes of structure factors from single crystals - $|F_o|$, σ_o . It is the most common case. Usually structure factors are reliable, uncertainties are not so.**
- 2) **Intensities or amplitudes are from “twinned” crystals**
- 3) **Amplitudes of structure factors are available for $|F^+|$ and $|F^-|$ - SAD case**
- 4) **Amplitudes of structure factors are available from multiple crystal forms**
- 5) **Complex structure factors are available – e.g. cryo-EM models**

Two components of target function

Crystallographic target functions have two components: one of them describes the fit of the model parameters into the experimental data and the second one describes chemical integrity (restraints).

Currently used restraints are: bond lengths, angles, chirals, planes, ncs if available, some torsion angles, reference structures

Crystallographic refinement

The function in crystallographic refinement has a form:

$$L(p) = wL_X(p) + L_G(p)$$

Where $L_X(p)$ is -loglikelihood and $L_G(p)$ is -log of prior probability distribution – restraints: bond lengths, angles etc.

It is one of many possible formulations. It uses Bayesian formulation. Other formulation is also possible.

-loglikelihood

-loglikelihood depends on assumptions about the experimental data, crystal contents and parameters. For example with assumptions that all observations are independent (e.g. no twinning), there is no anomalous scatterers and no phase information available, for acentric reflections it becomes:

$$L_X(p) = \frac{|F_o|^2 + |F_c|^2}{2s^2 + S} - \log(I_0(2|F_o||F_c|/(2s^2 + S)) + \log(2s^2 + S) + \text{const}$$

And for centric reflections:

$$L_X(p) = \frac{|F_o|^2 + |F_c|^2}{2(s^2 + S)} - \log(\cosh(|F_o||F_c|/(s^2 + S))) + 0.5\log(s^2 + S) + \text{const}$$


All parameters (scale, other overall and atomic) are inside $|F_c|$ and Σ

Note that these are loglikelihood of multiples of chi-squared distribution with degree of freedom 2 and 1


TWIN

The whole crystal: twin or polysynthetic twin?

twin

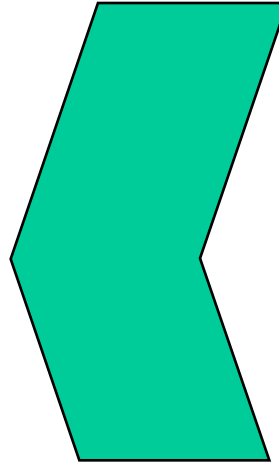

polysynthetic
twin

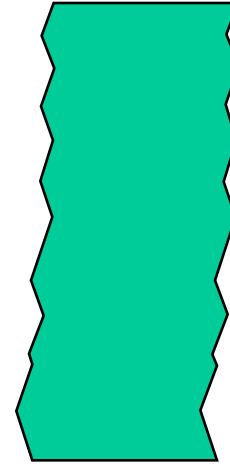


A single crystal can be cut
out of the twin:

yes

no

The shape of the crystal suggested that we dealt with polysynthetic OD-twin

Twin refinement: Group/subgroup

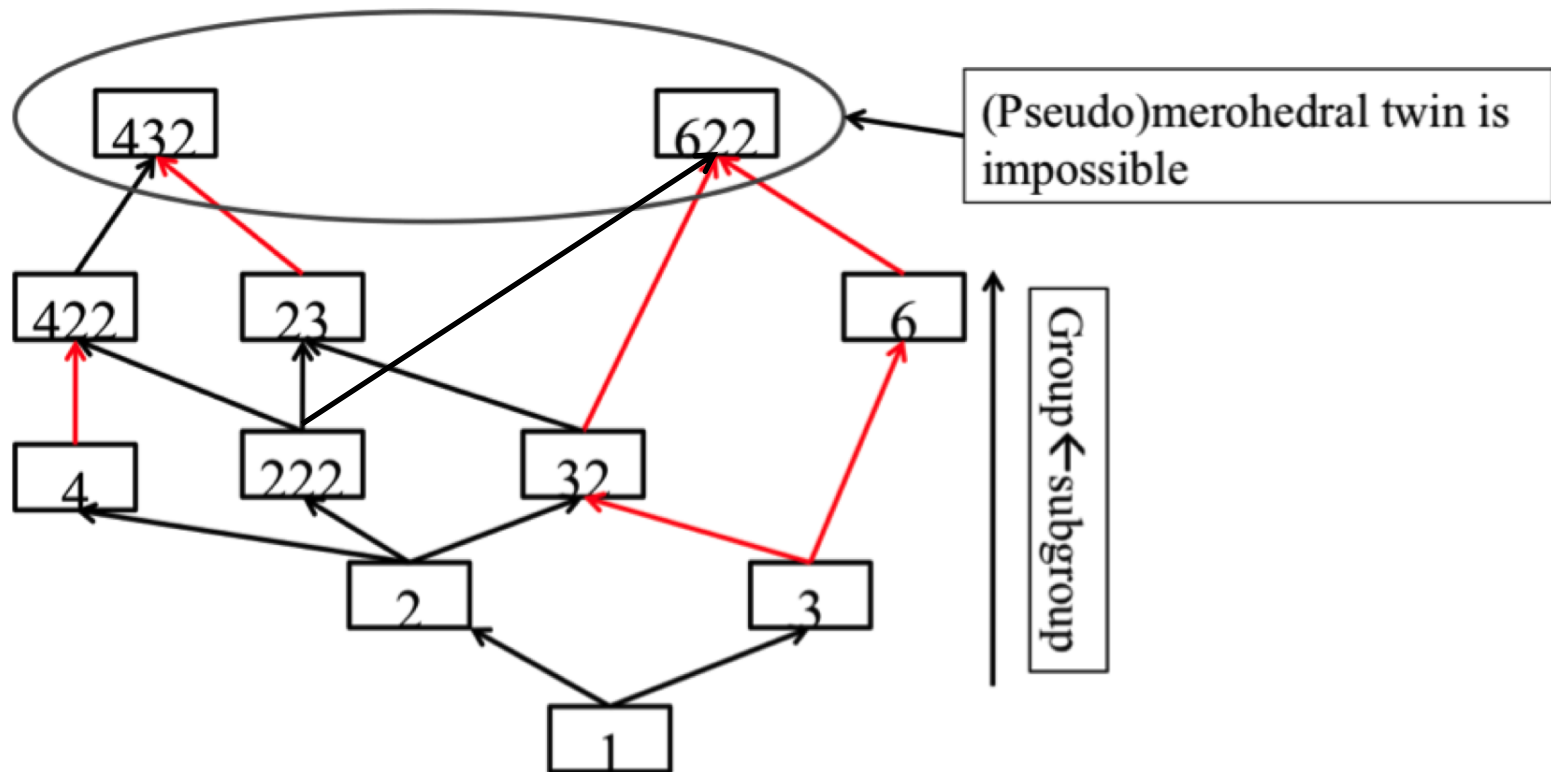

Red arrows: No constraints are needed, merohedral twin could happen
 Black arrow: Additional constraints on cell parameters are needed, pseudo merohedral twinning can happen

31

Twin refinement in REFMAC

Twin refinement in refmac (5.5 or later) is automatic.


- Identify “twin” operators
- Calculate “Rmerge” $(\sum |I_h - \langle I \rangle_{\text{twin}}| / \sum I_h)$ for each operator. If Rmerge > 0.44 remove this operator.
- Refine twin fractions. Keep only “significant” domains (default threshold is 7%): Twin plus symmetry operators should form a group

Intensities can be used

If phases are available they can be used


Maximum likelihood refinement is used

Where's the density for my ligand (2.15Å)?

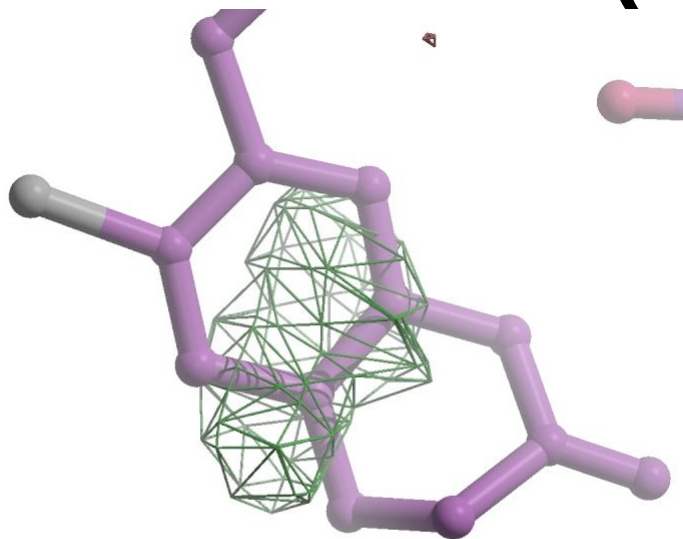

R-factor (R-free) 25.5% (26.9%) –
after initial rigid body and restrained
refinement.

Fo-Fc – 3 sigma

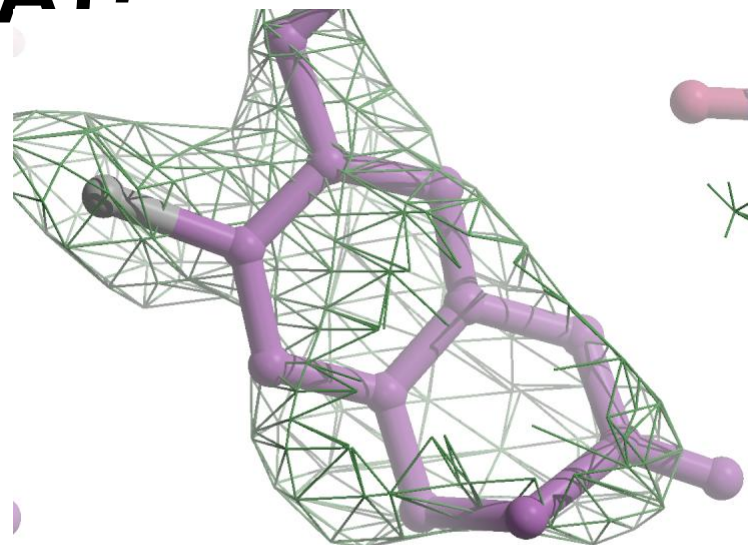

R-factor (R-free) 15.9% (16.3%) – re-
run restrained ref. with twin on
(refined twin fractions
0.6043/0.3957).

Fo-Fc – 3 sigma

Borrowed from B. Bax, GSK, Stevenage, UK

Overall R factor

Overall R-factor:

$$R_{work/free} = \frac{\sum_{work/free} \|F_o\| - \|F_c\|}{\sum_{work/free} \|F_o\|}$$

R-factor 1) depends on statistical property of data; 2) it is very sensitive to weights.

R-factor can be replaced by average R-factor. But it will still be dependent on statistical properties of data

Note: when sharpening is used then observed structure factors are multiplied by B-value dependent function: $\exp(-B s^2/4)$. After scaling F_c are also multiplied by the same function.


Twin: Few warnings about R values

Rvalues for random structures (no other peculiarities)

Twin	Modeled	Not modeled
Yes	0.41	0.49
No	0.52	0.58

Murshudov GN “Some properties of Crystallographic Reliability index – Rfactor:
Effect of Twinning” Applied and Computational Mathematics”, 2011:10;250-261

Rvalue for structures with different model errors:
Combination of real and modeled perfect twin fractions

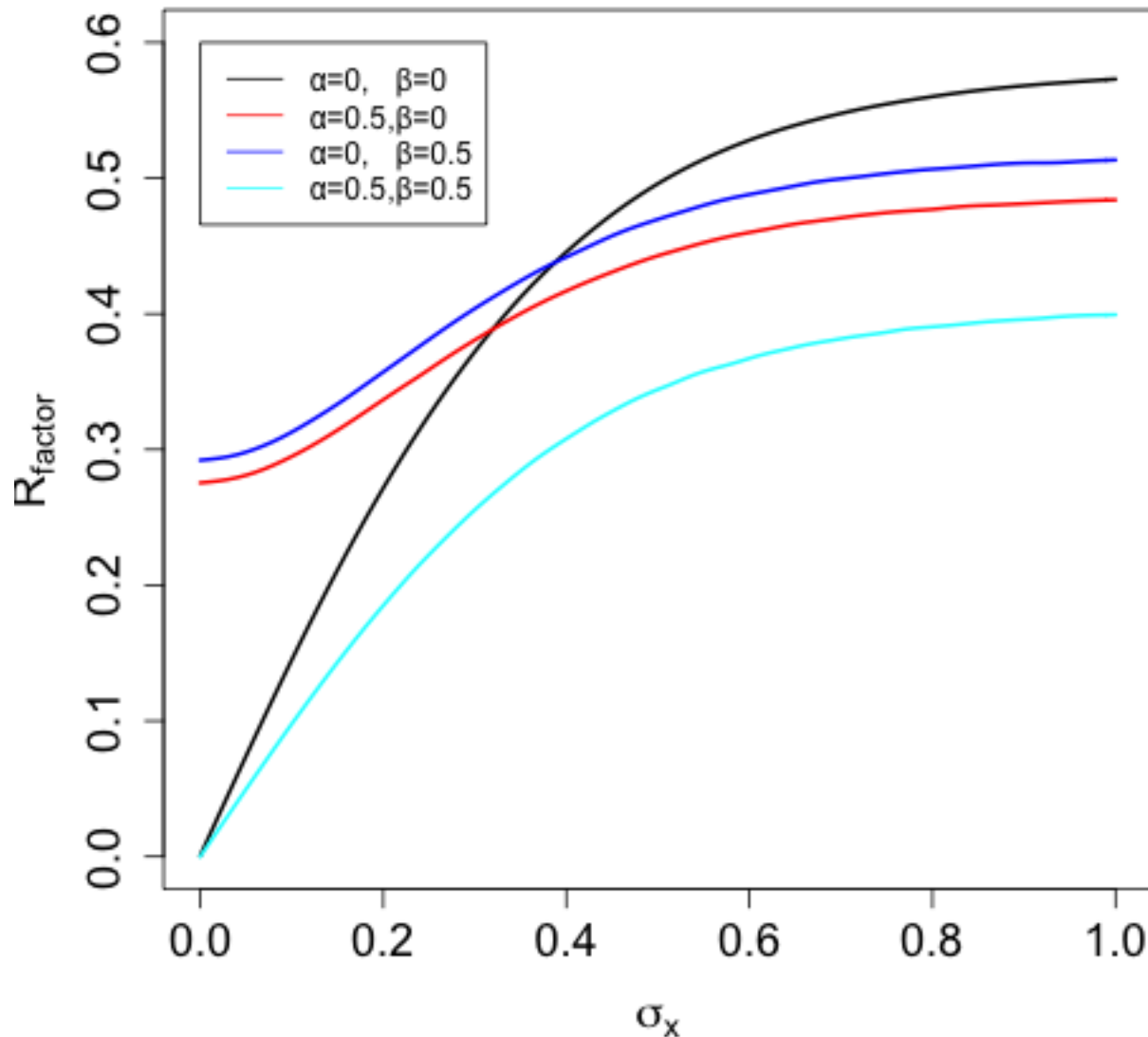

Overall R factor: average correlation

Average correlation coefficient seems to be less dependent on data properties and resolution dependent weights

$$CC_{average} = \frac{\sum_{i=1}^{reso:bins} N_i CC_i}{\sum_{i=1}^{reso:bins} N_i}$$

N_i number of reflections in resolution bin i

CC_i Correlation coefficient in resolution bin i

Example

Bvalue	R-factor	CC overall	CC average
0	26	92	80
50	25	94	80
100	24	95	80
-100	29	82	80

Parameters

Usual parameters (if programs allow it)

- 1) Positions x,y,z**
- 2) B values – isotropic or anisotropic**
- 3) Occupancy**

Derived parameters

- 4) Rigid body positional**
 - After molecular replacement
 - Isomorphous crystal (liganded, unliganded, different data)
- 5) Rigid body of B values – TLS**
 - Useful at the medium and final stages
 - At low resolution when full anisotropy is impossible
- 6) Torsion angles**

TLS

Rigid-body motion

General displacement of a rigid-body point can be described as a rotation along an axis passing through a fixed point together with a translation of that fixed point.


$$\underline{u} = \underline{t} + \underline{D}\underline{r}$$

for small librations

$$\underline{u} \approx \underline{t} + \underline{\omega} \times \underline{r}$$

$\underline{D}$ = rotation matrix

$\underline{\omega} \times$ = vector along the rotation axis of magnitude equal to the angle of rotation

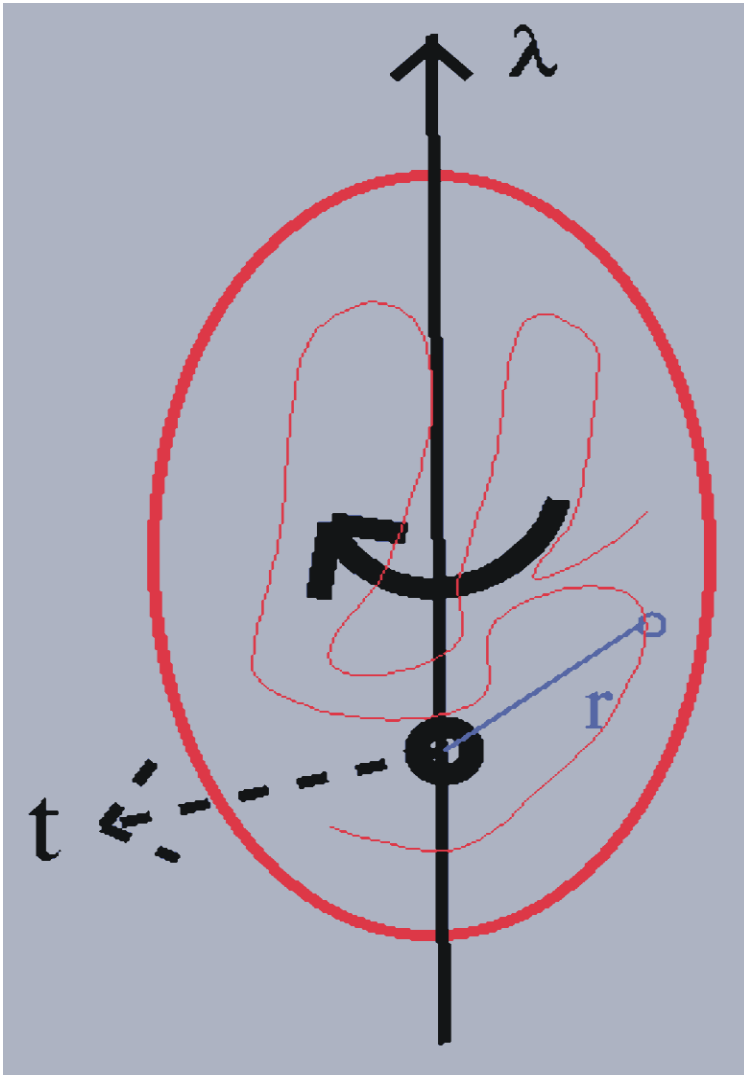

TLS groups

Rigid groups should be defined as TLS groups. As starting point they could be: subunits or domains.

If you use script then default rigid groups are subunits or segments if defined.

In ccp4i you should define rigid groups (in the next version default will be subunits).

Rigid group could be defined using TLSMD webserver:


<http://skuld.bmsc.washington.edu/~tlsmd/>

Bulk solvent

Method 1: Babinet's bulk solvent correction

At low resolution electron density is flat. Only difference between solvent and protein regions is that solvent has lower density than protein. If we would increase solvent just enough to make its density equal to that of protein then we would have flat density (constant). Fourier transformation of constant is zero (apart from F000). So contribution from solvent can be calculated using that of protein. And it means that total structure factor can be calculated using contribution from protein only

$$\begin{array}{ll} \rho_s + \rho_p = \rho_T & \Leftrightarrow F_s + F_p = F_T \\ \rho_s + k\rho_p = c & \Leftrightarrow F_s + kF_p = 0 \\ F_s = -kF_p & \Rightarrow F_T = F_p - kF_p = (1-k)F_p \end{array}$$

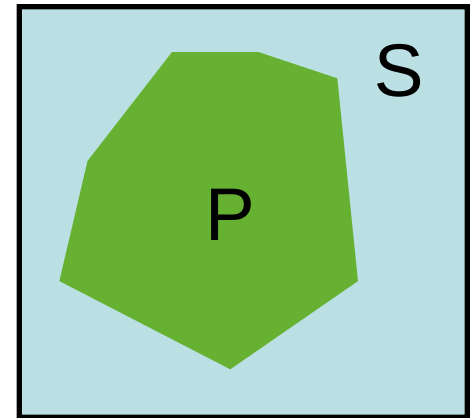

k is usually taken as $k_b \exp(-B_b s^2)$. k_b must be less than 1. k_b and B_b are adjustable parameters


Bulk solvent

Method 2: Mask based bulk solvent correction

Total structure factor is the sum of protein contribution and solvent contribution. Solvent region is flat. Protein contribution is calculated as usual. The region occupied by protein atoms is masked out. The remaining part of the cell is filled with constant values and corresponding structure factors are calculated. Finally total structure factor is calculated using

$$F_T = F_p + k_s F_s$$

k_s is adjustable parameter.

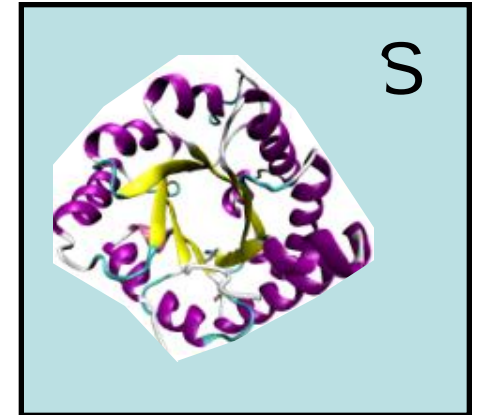

Mask based bulk solvent is a standard in all refinement programs. In refmac it is default.

Overall parameters: Scaling

There are several options for scaling:

- 1) Babinet's bulk solvent assumes that at low resolution solvent and protein contributors are very similar and only difference is overall density and B value. It has the form:

$$k_b = 1 - k_b \exp(-B_b s^2/4)$$

- 1) Mask bulk solvent: Part of the asymmetric unit not occupied by atoms are assigned constant value and Fourier transformation from this part is calculated. Then this contribution is added with scale value to "protein" structure factors. Total structure factor has a form:

$$F_{\text{tot}} = F_p + s_s \exp(-B_s s^2/4) F_s.$$

- 1) The final total structure factor that is scaled has a form:

$$s_{\text{aniso}} s_{\text{protein}} k_b F_{\text{tot}}$$

Map calculation

- After refinement programs give coefficients for two type of maps: 1) 2Fo-Fc type maps. 2) Fo-Fc type of maps. Both maps should be inspected and model should be corrected if necessary.
- Refmac gives coefficients:

$2 m F_o - D F_c$ – to represent contents of the crystal

$m F_o - D F_c$ - to represent differences

m is the figure of merit (reliability) of the phase of the current reflection and D is related to model error. m depends on each reflection and D depends on resolution. Unobserved reflections are replaced by DF_c .

If phase information is available then map coefficients correspond to the combined phases.

Available refinement programs

- SHELXL
- CNS
- REFMAC5
- TNT
- BUSTER/TNT
- Phenix.refine
- RESTRAINT
- MOPRO
- XD
- MAIN

What can REFMAC do?

- **Simple maximum likelihood restrained refinement**
- **Twin refinement**
- **Phased refinement (with Hendrickson-Lattmann coefficients)**
- **SAD/SIRAS refinement**
- **Structure idealisation**
- **Library for more than 8000 ligands (from the next version)**
- **Covalent links between ligands and ligand-protein**
- **Rigid body refinement**
- **NCS local, restraints to external structures**
- **Helical, point group NCS constraints**
- **TLS refinement**
- **Fit into EM map**
- **Map sharpening**
- **etc**

What and when

- Rigid body: At early stages - after molecular replacement or when refining against data from isomorphous crystals
- “Jelly” body – At early stages and may be at low resolution
- TLS - at medium and end stages of refinement at resolutions up to 1.7-1.6Å (roughly)
- Anisotropic - At higher resolution towards the end of refinement
- Adding hydrogens - Higher than 2Å but they could be added always
- Phased refinement - at early and medium stages of refinement
- SAD – at the early stages
- Twin – when you are sure that crystal is twinned
- NCS local – always?
- Ligands - as soon as you see them
- What else?

Conclusion

- Information about ligands and their chemistry should be used in refinement
- NCS restraints are useful tool
- External restraints can be a powerful tool for reliable atomic model derivation.
- “Jelly” body can be very powerful at the early stages of refinement
- Twin refinement improves statistics and occasionally electron density: Rfactors may be misleading
- Refinement is just one step in X-ray structure analysis – it is often used as part of phase improvement and model building

Acknowledgment

Alexei Vagin

Andrey Lebedev

Rob Nicholls

Fei Long

Paul Emsley

David Tourigny

Andrea Thorn

Pavol Skubak

Raj Pannu

Lamzin group, Hamburg

Perrakis group, Amsterdam

CCP4, LMB people

REFMAC is available from CCP4 and from our website:

<http://www2.mrc-lmb.cam.ac.uk/groups/murshudov/>